- Never use white dross.
- This manual ONLY uses attachments recommended by the manufacturer.
- Use this product only as described or as assembled in this manual, or other children's related products.
- Keep your child away from the product when it is running.
- Drop the cord away from the plug - NOT THE CORD.
- Check if grasping the plug - NOT THE CORD.
- DO NOT connect any hand or pull the cord (pulling).
- Refer to the Advanced Service Center for repairs.
- If it is not working properly, if it has been dropped.
- Never operate the unit if it has a damaged plug.

DANGER!

To reduce the risk of burns, electrocution, fire, or injury:

- ALWAYS unplug this product immediately after using.
 - DO NOT place or store the unit where it can fall or be pulled into a tub, sink, or toilet.
 - DO NOT reach for a unit that has fallen into water.
 - Unplug immediately. Do not reach into the water under any circumstances. The unit has been designed for use in a 20 volt AC.
 - The unit is designed for use in a 20 volt AC.
 - DO NOT submerge the unit in water or any other liquid.

.....

The Interplak difference

The Interplak® Dental Water Jet was specifically designed to help you maintain a high level of oral hygiene and improved oral health quickly and easily at home. It flushes food debris from those hard-to-reach places around implants, orthodontic braces and crowns. It is gentle and effective on all orthodontic, crown, and bridge work.

FCC Information to the user

This device complies with Part 18 of the FCC Rules. Operation is subject to the following two conditions:
(1) this device may not cause harmful interference, and (2) this device must accept any interference that may cause undesired operation.

Emissions Note: This equipment has been tested and found to comply with the limits pursuant to Part 18 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and if not installed and used in accordance with instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.


If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

- ### ***Cleaning and Care***
- If, in the course of use, the valve (8) on the water tank begins to leak, this can be taken care of by cleaning the valve seat. To do this, push out the valve body using the nozzle tip, clean the valve seat and replace the valve.

Cleaning and Care

If, in the course of use, the valve (8) on the water tank begins to leak, this can be taken care of by cleaning the valve seat. To do this, push out the valve body using the nozzle tip, clean the valve seat and replace the valve.

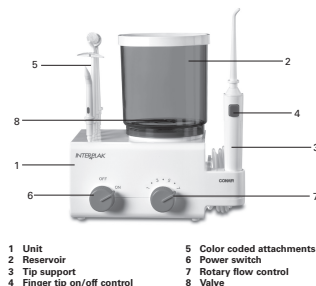© 2009 Cengage Learning. All Rights Reserved. May not be copied, scanned, or duplicated, in whole or in part. WCN 02-200-203

Instructions For Use

Remove the water reservoir (2) and fill with lukewarm water. Reinstall the reservoir firmly onto the unit (1). Prior to initial use, the unit should be rinsed by letting a full reservoir of water pass through it. If the unit is to be used with a mouthwash solution, it should be rinsed with clean water after each use to remove any traces of mouthwash. If the unit has not been used for a long period of time, for reasons of hygiene it is advisable to follow the directions for initial use, letting a full reservoir of water pass through. To remove trapped air upon initial operation, turn the rotary flow control (7) all the way to the right and then to the left, three times.

Prior to filling the water reservoir, make sure water tip is firmly placed on and in the stand. The Water Jet is supplied with 5 attachments. Remove the tip support (3) and install selected tip:

Standard Tip (2pcs): removes trapped particles from between teeth.
Subgingival Tip (1pc): uses a soft rubber tip to massage the gums and safely deliver solution into the gingival pockets below the gumline.
Brush Head Tip (1pc): loosens debris and larger food particles before flushing them away – ideal for all dental work including fixed and removable appliances.


Tongue Cleaner (1pc): place the tongue cleaner at farthest end of tongue and move forward. Repeat to cover entire tongue. After switching on (6) the unit, the rotary flow control (7) can be used to regulate the water pressure. You should begin with low pressure and gradually increase it to the desired level. Hold the fingertip on/off control (water-stop button) (4) down and, after switching on the unit, place the tip in your mouth and lean over the sink. When you release the water-stop button, the cleaning and massage process is begun. The direction of the water jet in your mouth can be controlled by twisting the tip to the desired position on the tongue. The water jet is directed away from the tongue so that the water can run back into the sink. Direct the water jet toward your teeth, especially between your teeth, in order to rinse away stubborn remaining food. The massage you receive while cleaning your tongue also massages your gums. The water jet is directed toward your gums and the water pressure is increased to remove debris, relax and soothe your gums.

Welcome to a New Level of At-Home Dental Care

The Interplak® Water Jet enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

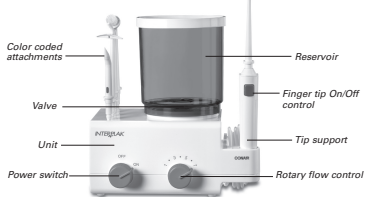Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the Water Jet can most effectively meet your specific oral needs.

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