ALARME DE LIT SANS FIL

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WIRELESS BED ALARM

Thank you for purchasing an Ideal Security product.

The 6311 Wireless Bed Alarm includes a weight sensing pad (sensor), a wireless transmitter, and a wireless receiver alert.

We take customer service seriously and pride ourselves on our relationships with our customers. If you encounter any difficulties during the set-up or have any questions, please see the customer service section of this document.

FEATURES

WIRELESS RECEIVER ALERT FEATURES

- Six different alarm tones, for distinguishing which pad was triggered if multiple pads are set up
- Silent (lights only) alarm mode available
- Power: 3 x C Alkaline batteries required (not included)
- Optional AC adapter (not included) available for purchase
- Battery life in standby mode: 8 months approx.
- Battery-low LED indicators
- Sound level (at 1 ft / 30 cm): 95 dB
- Signal transmission range: 131 ft (40 meter), line of sight
- Wireless signal frequency: 433.92 MHz
- Mountable using the keyhole slots on the rear. Hardware not included
- Size: 5" x 3-3/4" x 1-1/2" (128 x 95 x 38 mm)
- Indoor use only
- Wireless receiver alert can link up to six weight sensing pad with wireless transmitters

ALARM FEATURES

By default, the transmitter and receiver are factory-linked to the #5 sound (see table in the next section).

Up to 6 transmitters can be linked to a single receiver. Each transmitter will be wirelessly connected to one of the six sounds on the receiver, having its own unique sound. Each sound can be linked to one transmitter only. Once a transmitter is linked to a sound, the sound is then removed from the available sounds.

The receiver also has a silent mode available, which flashes light with no sound. This alert will last as long as the selected sound on the receiver. (i.e. if the receiver has sound #4 selected, the light will flash for 30 seconds).

A single transmitter can also be linked to multiple receivers, triggering alerts on all receivers at the same time.

Please note that multiple receivers should be placed at least 3 ft (1 m) apart from each other.

Listed below are the six different sounds in order:

Sound	Duration
1. Ding-dong	2 cycles
2. Di-di	10 seconds
3. Di-di-di	30 seconds
4. Di-du	30 seconds
5. DididiDududu	Forever, batteries will last up to 36 hours
6. Wawa siren	Forever, batteries will last up to 36 hours

Silent mode - lights only. Same as selected sound (See SOUND OFF MODE in OPERATION SECTION).

WIRELESS TRANSMITTER FEATURES


- Power: 2 x AAA Alkaline batteries required (not included)
- Replaceable batteries: Unlike other mats, the batteries in your 6311 transmitter are easily replaceable
- Battery life in standby mode: 1 year approx.
- Battery low indicator
- Transmission distance: Up to 40 m (130 ft), line of sight
- Wall mountable with the included double-sided tape or using screws (not included)
- Long lead wire with connector: 10 ft for easy placement
- Size: 2-7/8" x 1-1/4" x 1" (72 x 31 x 25 mm)
- Indoor use only
- FCC approved device

WEIGHT SENSING PAD FEATURES


- Size: 30" x 14" x 3/16" thick (762 x 355 x 4.75 mm)
- Sensitivity: 1 psi
- IP44 water splash resistant pad
- Indoor use only
- Easy to replace (replacements available for order)
- Long lead wire with connector: 6 ft

DESCRIPTION OF COMPONENTS

RECEIVER

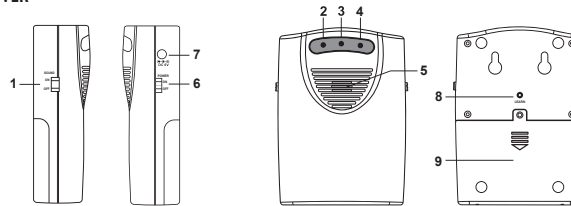
- | | |
|--|--------------------------|
| (1) Sound on/off switch | (5) Buzzer |
| (2) Transmitter battery low/trigger indication LED | (6) Power on/off switch |
| (3) In-use indication LED | (7) DC jack |
| (4) Receiver battery low/trigger indication LED | (8) Learn ID code button |
| | (9) Battery cover |

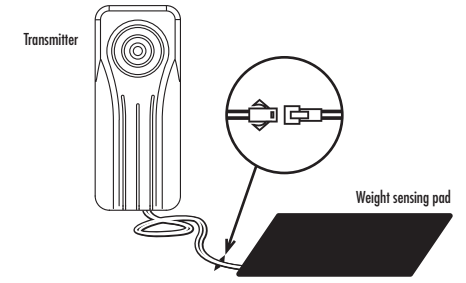POSITIONING


We recommend placing the weight sensing pad under the shoulders or hips, where the most consistent contact is made with the bed. Generally, the shoulders provide the most consistent performance, but you may find the pad more comfortable if placed under the hips. If you toss and turn frequently at night, you may experience more false alarms. Try different positions to optimize for comfort and reliability.

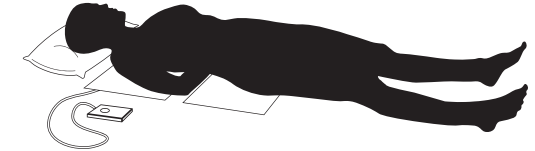For safety, we recommend hiding away the wires under the bed structure and/or taping the wires onto a wall or floor. A loose wire may present a fall hazard.

SET-UP AND INSTALLATION

RECEIVER

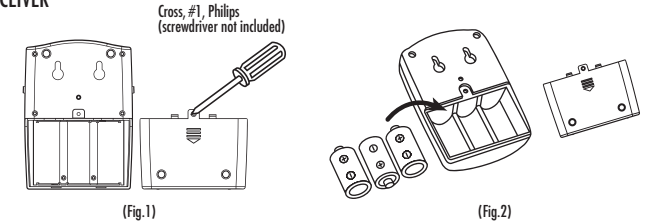
1. Remove screw from back of the receiver cover and open the panel. (Fig.1)
2. Install three "C" batteries in the receiver. (Fig.2)
3. You can connect an AC adapter (6V DC 500mA (not included); plug length: 12mm; inside positive ID: 2.1mm; outside negative OD: 5.5mm) to the alert, and install the batteries as backup in case of power failure.

WIRELESS TRANSMITTER

1. Using #1 screwdriver (not included), remove the screw from the bottom of the transmitter and remove cover.
2. Install two AAA batteries (not included).

Note: There is no ON/OFF switch. If you want to turn the alarm off you can do so on the receiver.

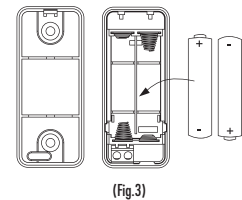OPERATION

SET-UP YOUR WIRELESS BED ALARM QUICKLY WITH THE FACTORY LINKED SOUND

The shipped transmitter and receiver were factory linked to sound #5 already.

1. Plug the connectors of the transmitter and weight sensing pad together.
2. Place the weight sensing pad on a big flat surface that can support the entire pad.
3. Place the transmitter away from the pad, not on the pad. You can position the transmitter using the included double-sided tape or you can simply place it on a surface close by (i.e., bedside table).
4. On receiver, slide the POWER and SOUND switches to the ON position.
 - You will immediately hear a beep. Now the system is in standby/ready mode.
 - The ARMED middle LED light will flash every 3 seconds indicating the system is ready.
5. The transmitter is armed when the pad senses weight. To test, sit on pad.
6. A signal will trigger when the weight is removed from pad. To test, sit up and off of pad.

Sound ON mode: When the sensor is triggered, the alarm will sound, and lights will flash.

Sound OFF mode: When the sensor is triggered, only lights will flash.

To stop the alarm: If sounds 1 - 4 are selected, the alarm will end on its own. If a "forever" alarm is selected (#5 or #6), slide the POWER switch to OFF position.

To restart the alarm: If sounds 1 - 4 are selected, the alarm is ready on its own. If a "forever" sound is selected (#5 or #6), slide the POWER to the OFF position then to the ON position.

BATTERY LOW INDICATOR

When the batteries in the receiver alert are low, the right LED (#4 description of components receiver) will flash once every four seconds.

When the batteries in the transmitter are low, the integrated LED will flash slowly, 3 seconds ON, one second OFF. Also, the #2 LED in the receiver will flash.


Change the batteries immediately when battery is low.

LINKING A TRANSMITTER AND RECEIVER

The shipped transmitter and receiver were factory linked to sound #5 already.

These instructions are for linking a new transmitter to a receiver.

1. Plug the connectors of the transmitter and weight sensing pad together.
2. Place the weight sensing pad on a big flat surface that can support the entire pad.
3. Place the transmitter away from the pad, not on the pad, you can position the transmitter using the included double-sided tape or you can simply place it on a surface close by (i.e., bedside table).
4. On receiver, slide the POWER switch to ON position.
5. On receiver, slide the SOUND switch to ON position.
6. Sit on the pad. You will sit on the pad during the whole process unless otherwise instructed.
7. On receiver, slide the POWER switch to OFF position.
8. Sit on the pad continuously for at least 15 seconds.
9. On receiver, slide the POWER switch to ON position.
10. Locate the LEARN button on the back of the receiver, push it for 1 second then release (Fig. 4).

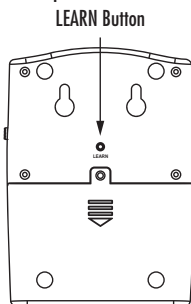(Fig.4)

11. You will hear a sound. If you like the sound, go to next step, otherwise, push the button again for 1 second then release to hear the next sound, you can cycle this way through all available sounds.
12. Once you hear the sound you like, push the learn button continuously for 3 seconds until you hear a beep. Release the button. Now the receiver is ready to accept the new transmitter.
13. Stand up from the pad.
14. You will hear one long beep sound. That means the linking (learn) is successful.
15. If you hear 4 short beeps, that means the linking has failed. Slide POWER to OFF on the receiver then restart the procedure from step 1.
16. Test the system again.
17. On the receiver, make sure POWER and SOUND are set to ON.
18. Sit back on the pad for 5 seconds.
19. Stand up and off the pad.
20. You will hear the sound you selected for the duration associated with the sound. That means the system works.
21. For sound #5 and #6 where the sound continues indefinitely, you can turn POWER to OFF on the receiver and turn it back ON. This will not erase the linked transmitter.

Now you have successfully linked a new transmitter to a receiver.

ASSIGNING A NEW SOUND OR REMOVING A TRANSMITTER

Follow the instructions below if you want to:

- Change the default sound #5, you can choose a new sound and link to the transmitter
- Unlink (remove) a transmitter
- Relink a transmitter to receiver after system has been inactive for a period of time

1. On receiver, make sure POWER is set to OFF.
2. Do step 1 to 9 as described in the previous section.
3. Locate the LEARN button on the back of the receiver, push it continuously for 12 seconds until you hear two long beeps, release the button. This will erase all linked transmitters from the receiver. You cannot erase only one transmitter.
4. Do step 8 to 21 in the previous section.

CLEANING, CARE AND MAINTENANCE

While the pad is splash resistant, it is not water-proof and is not designed for submersion in water.

We recommend placing the weight sensing pad under a washable bedding or incontinence pad. Be sure to test the function before using. If the protective layer is too heavy, it may prevent the weight sensing pad from properly detecting weight changes.

To clean the pad, we recommend using disinfectant or anti-bacterial wipes.

If soiled, you may also clean the pad by using a mild soap in a spray bottle. Lightly spray the pad on both sides, taking care not to spray near the connector jack or where the wire enters the pad. Wipe down with a dry sponge or towel.

DO NOT submerge the pad in water or spray directly with water from a shower, hose, or faucet.

DO NOT use the pad when folded, as applying pressure could damage the internal circuitry.

Take the batteries out when the device will be out of use for a long period.

IMPORTANT NOTES

- Before installing permanently, test the reception range of the receiver. If the receiver does not pick up the signal from a transmitter, try another location closer to the transmitter.
- The components are designed for indoor use only. Do not install the receiver or transmitter where they will be exposed to direct sunlight, rain, or extreme cold environment.
- The sound coming from the receiver can be very loud. Never put the speaker close to your ear.
- For the most up to date installation instructions, videos and troubleshooting information consult our web site, www.idealinc.com
- This wireless bed alarm can provide additional security for your family and loved ones if utilized properly. However, its use cannot guarantee complete protection against falls or injury. Ideal Security cannot assume any responsibility for any damages as a result of installing this or any other Ideal Security product.

WARRANTY

This unit is warranted to be free of defective materials and workmanship for a period of one year from date of purchase. Please keep your bill of sale. In case of malfunction, this unit will be replaced at no charge by Ideal Security Inc. Warranty begins on date of consumer purchase. Replacement as described above shall be the sole remedy to the consumer and there shall be no liability on the part of Ideal Security Inc. for any special, indirect, incidental or consequential damages, including but not limited to any loss of business or profits whether or not foreseeable. For more information, visit www.idealinc.com/policies

CUSTOMER SERVICE

Thank you for your purchase.

We are dedicated to customer service and will work to resolve any difficulties to ensure customer satisfaction.

If you have any questions about the setup or function of this product, please contact us before returning it.

Tel.: 1-800-361-2236

Email: info@idealinc.com

Support and Contact: www.idealinc.com/help

FAQ and Videos: help.idealinc.com

IDEAL SECURITY INC.

LASALLE, QC • PRESCOTT, ON

www.idealinc.com

This device complies with Part 15 of the FCC Rules and RSS-210 of industry Canada.

1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesired operation of the device.

Warning: Changes or modifications made to this equipment not expressly approved by the party responsible for compliance may void the user's authority to operate this equipment.