

VitalStim® Therapy

Can Help Patients With Dysphagia

Dysphagia: a debilitating condition.


Dysphagia is defined as difficulty with swallowing. It can occur as a result of a stroke or other neurological disease, normal aging or after a long period of inactivity. An estimated > 16 million US, 30 million European and 10 million Japanese elderly citizens suffer from dysphagia.¹

People who don't eat or drink enough are more likely to become ill or experience longer recovery times. In addition, patients with dysphagia carry a large psychological burden; many develop depression and withdraw from society because sharing of food and drink is such a big part of daily life.

The management of this condition consists mainly of diet modification to decrease choking episodes; a feeding tube may be inserted in the stomach.

In the last few years, new treatment options have become available – especially the use of electrical stimulation. This exciting treatment tool, VitalStim® Therapy, is showing good outcomes in many patients. It is simple for trained clinicians to administer and is pain free for the patient.

1. Clavé P, Shaker R. Dysphagia: current reality and scope of the problem. Nat Rev Gastroenterol Hepatol. 2015 May; 12(5): 259-70

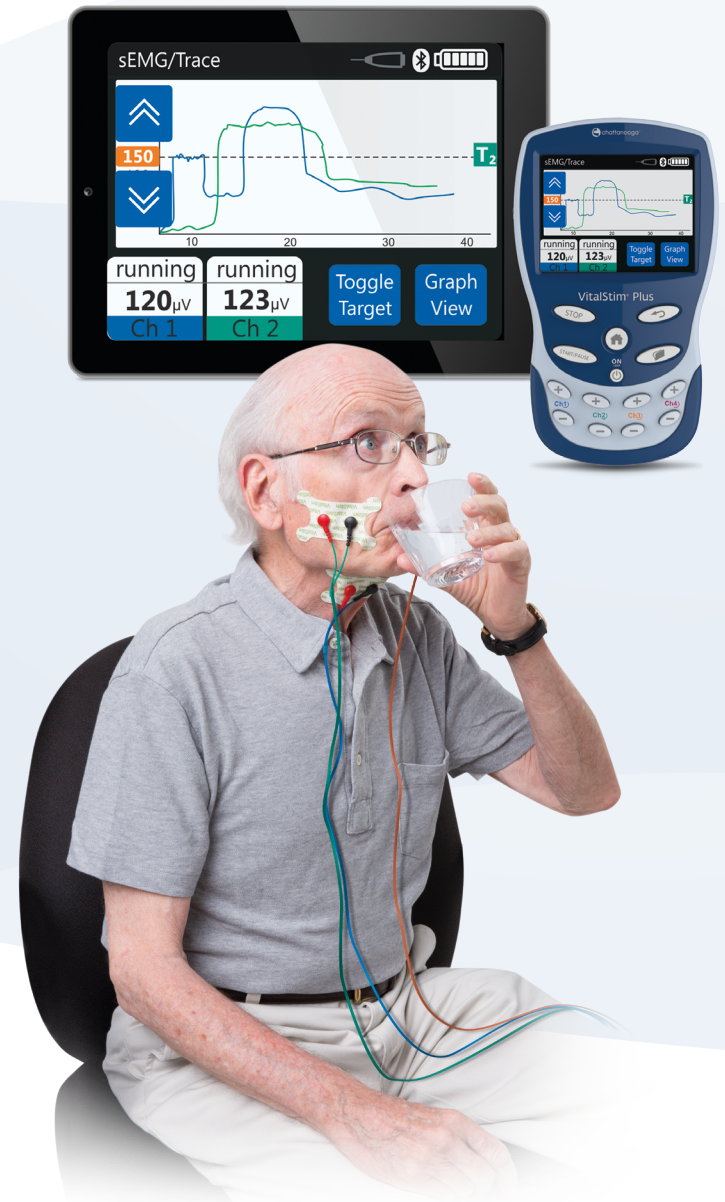

VitalStim® Plus

Electrotherapy and sEMG Biofeedback System



When you can't swallow



 **MOTION IS MEDICINE**⁺
DJO Global | 1a Guildford Business Park | Guildford
Surrey | GU2 8XG | UNITED KINGDOM
www.chattgroup.eu

© 2016 DJO - 000 - 00004195-EN - Rev A

Questions & Answers about
Treating Dysphagia with Electrical Stimulation



VitalStim® Therapy:

Treating Dysphagia with Electrical Stimulation

What To Expect When You Are Receiving VitalStim Therapy

If you've been diagnosed with dysphagia, you probably have a lot of questions. Dysphagia is defined as difficulty in swallowing. Patients with dysphagia may exhibit signs of penetration or aspiration (food or liquid getting into the lungs), muscle weakness and difficulty managing a normal diet. An impaired swallow often results in feeding restrictions, which in turn contribute to a decline in quality of life.

Who Is A Typical Patient?

Typical patient categories include, but are not limited to:

- Stroke
- Various other neuromuscular disease processes (e.g., Parkinson's, ALS, etc.)
- Head and neck cancer (post radiation) and/or surgery
- Deconditioning as a result of age or disease

Am I A Candidate For VitalStim Therapy?

If you show signs of aspiration or have difficulty managing your diet, you may be a candidate for dysphagia therapy. Look for one or more of the following signs and symptoms:

- Coughing/clearing throat after swallowing
- Abnormal volitional cough
- Decreased voice quality (wet, hoarse, weak)
- Recurring chest infections
- Multiple swallows or special maneuvers required to clear throat
- Difficulty completing a meal
- Feeling of food being stuck in the throat
- Modified diet required (thickening, pureeing food; soft solids)
- Difficulty initiating a swallow
- Spilling food or liquid from lips and/or drooling



Are There Any Contra-Indications For VitalStim Therapy?

There are few contraindications for receiving this therapy. With a complete medical history, your dysphagia therapist and medical team will make the decision with you about whether or not to use this specific therapy intervention.

Who Provides This Therapy?

VitalStim Therapy is performed by a dysphagia therapist. This is usually a speech language pathologist or physical therapist, or this can be your physician. In some countries a prescription by your physician may be required to receive VitalStim Therapy. The therapist will perform a clinical evaluation prior starting the therapy. At this time, information about your past and present medical history will be taken. Information about your current cognitive status, functional swallowing ability and your own perception will be recorded. Usually, at this time, other exams will be scheduled either with x-ray, ENT (ear nose throat) or GI (gastrointestinal) as appropriate. The dysphagia therapist will decide whether VitalStim Therapy is indicated after the initial evaluation is concluded.

How Does VitalStim Therapy Work?

VitalStim Therapy is designed to help swallowing through the use of Neuromuscular Electrical Stimulation (NMES). Electrical stimulation is used to aid muscle strengthening to rehabilitate the swallow. The sensory stimulation that is provided also helps muscle recruitment and swallowing function.

A VitalStim Therapy program is typically administered 3 to 5 days per week for 30-60 minutes. Patients can expect their treatment program to last a few weeks, sometimes longer.

A typical treatment session begins with the dysphagia therapist cleansing the skin with a special alcohol swab. Sometimes, moist heat and massage will be applied to improve the effects of the treatment. The dysphagia therapist will then attach the electrodes to the neck. Once turned on, the VitalStim device provides stimulation to the target muscles. At this time, your therapist will coach you through active swallowing exercises that include: "Swallow hard, swallow fast" and "swallow in a single contraction". These swallows will include food and liquids that are deemed appropriate by your dysphagia therapist.



VitalStim® Plus

Electrotherapy and sEMG Biofeedback System

Preventing Complications:

- Continue working with your physician to manage your reflux with proper medications and lifestyle changes.
- Follow all diet modifications as instructed by your dysphagia therapist.
- Report any changes in your swallow ability to your dysphagia therapist (i.e., increased coughing, choking, etc).
- Follow your home program.
- Track your meals and oral intake at home in a food log as instructed by your dysphagia therapist.
- Report any changes in your respiratory status to your dysphagia therapist (i.e., fever, cough, chest congestion, etc.).
- Attend all scheduled therapy sessions.
- Work with your therapist to modify your schedule if conflicts arise.

To learn more about VitalStim Therapy or to find an official Chattanooga VitalStim Plus Therapist near you, visit www.chattgroup.eu/vitalstim