

WE HOPE TO SEE YOU THERE!

A Systematic Review of Needleless Connector Function and Occlusion Outcomes

Monday, Nov. 10th, 2025

6:30 PM – 8:30 PM

Cooper's Hawk,
Virginia Beach

Presenter: Patrick Kiley, RN BSN



THE SOLUTION: NEXUS TKO® (TO KEEP OPEN) ANTI-REFLUX TECHNOLOGY


NEXUS TKO® REQUIRES NO CHANGE IN PRACTICE OR CLAMPING SEQUENCE AND PROVIDES BLOOD REFLUX PROTECTION 24/7™.

NEXUS TKO®-6P



NEEDLELESS CONNECTOR

ANTI-REFLEX VALVE

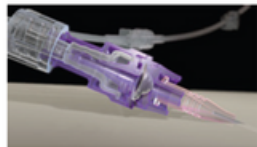

Cross sectional view of Nexus®TKO

HOW NEXUS TKO® WORKS

The Nexus TKO® anti-reflux valve opens and closes based upon fluid pressure. A micro slit in the valve opens towards the patient with fluid pressure from sources such as an IV or syringe pump, IV bag, or medication syringe. When the fluid pressure drops the Nexus TKO® valve instantly and automatically closes preventing blood reflux into the catheter.