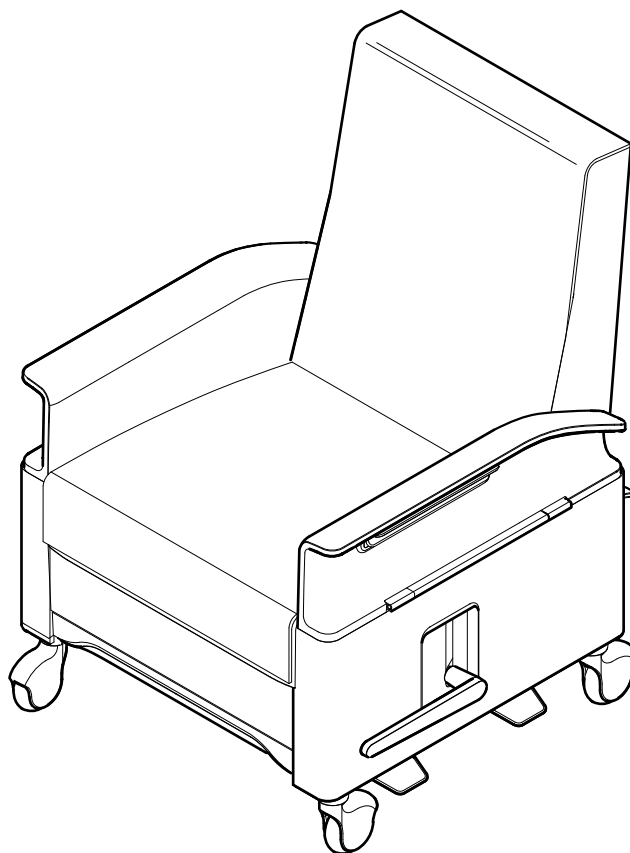


If you have a problem, question, or request, call
your local dealer, or Steelcase Line 1 at
888.STEELCASE (888.783.3522)
for immediate action by people who want to help you.


(Outside the U.S.A., Canada, Mexico, Puerto Rico,
and the U.S. Virgin Islands, call: 1.616.247.2500)
Or visit our website: www.steelcase.com

© 2012 Steelcase Inc.
Grand Rapids, MI 49501
U.S.A.

Printed in U.S.A.


Remove Pull Bar

1. Remove (2) set screws.
2. Pull outward to remove pull bar.
3. Loosen and remove bolts, washers, and extensions.

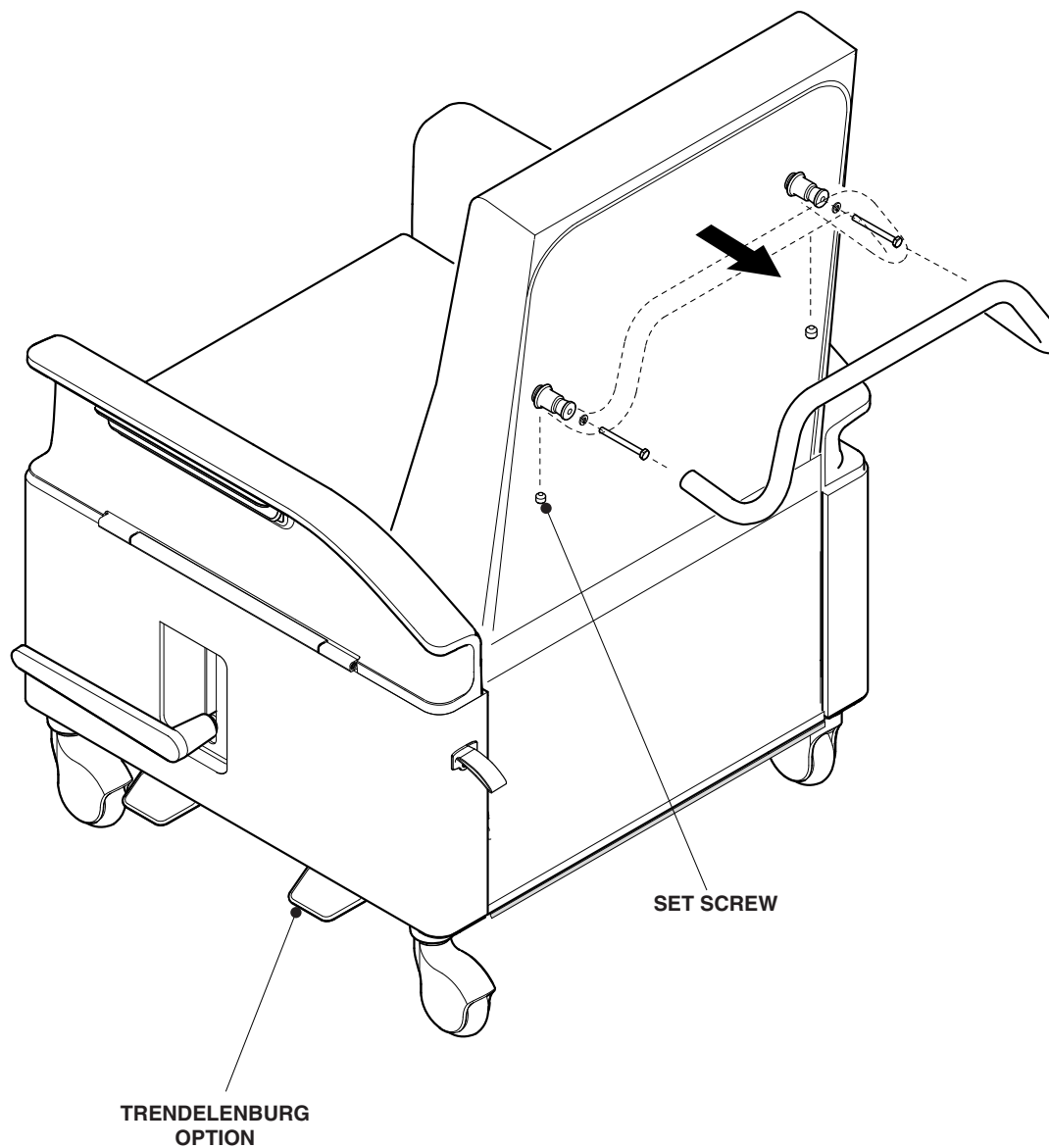


Bins for components removed



Steel

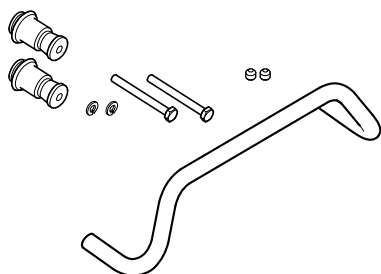



Remove Seat Back

- 1. Detach Velcro
- 2. Remove bayonet screws.
- 3. Pull upwards to remove seat back.
- 4. Cut off upholstery from seat back.
- 5. Pull foam off back frame (next page).

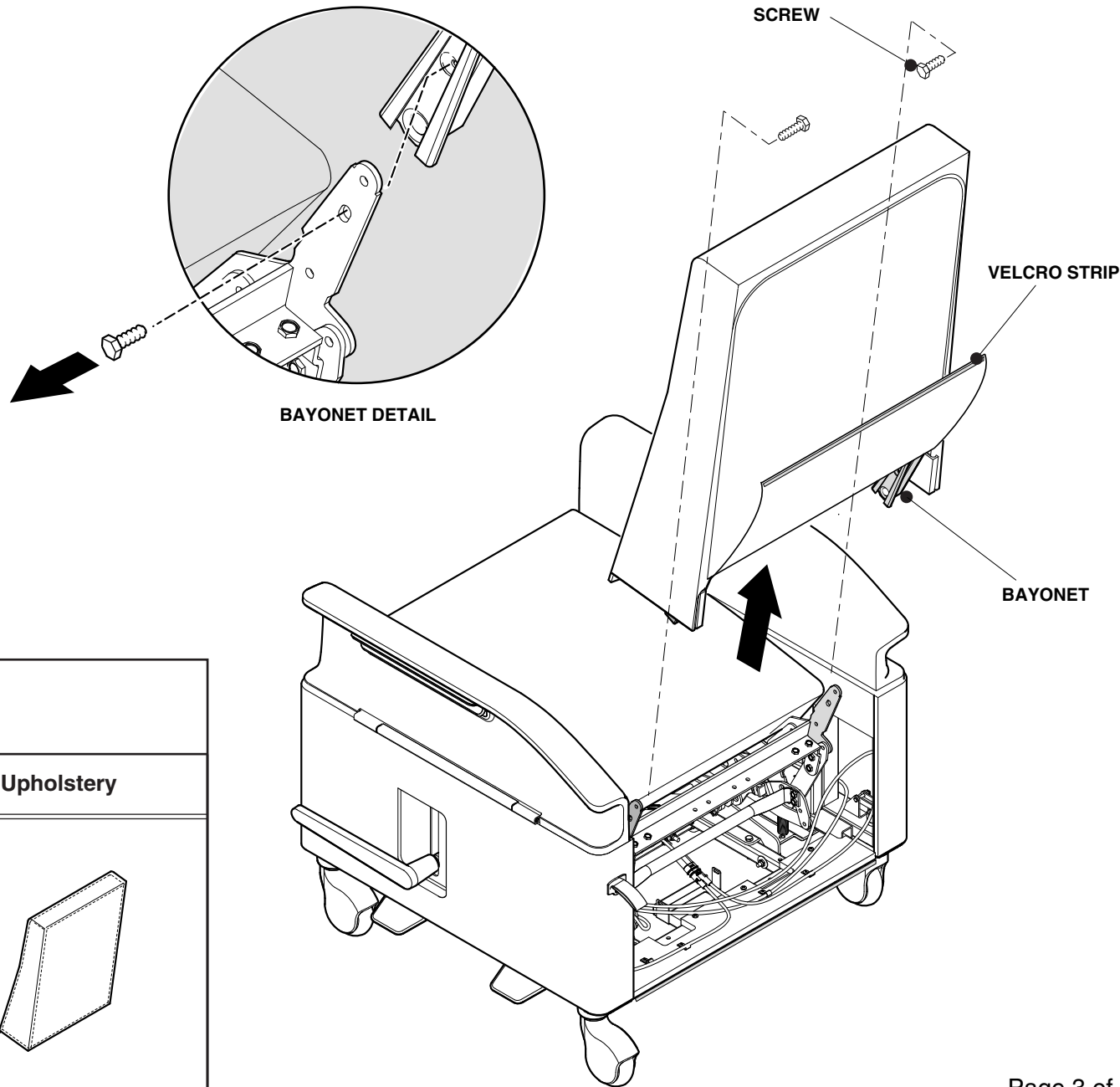

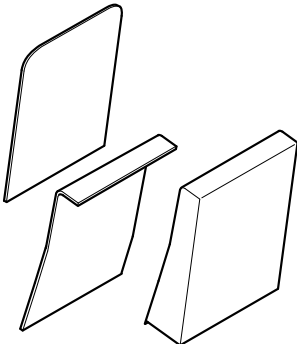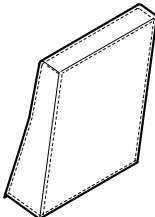Bins for components removed		
Steel	Foam	Upholstery
		


Disassemble Back Support

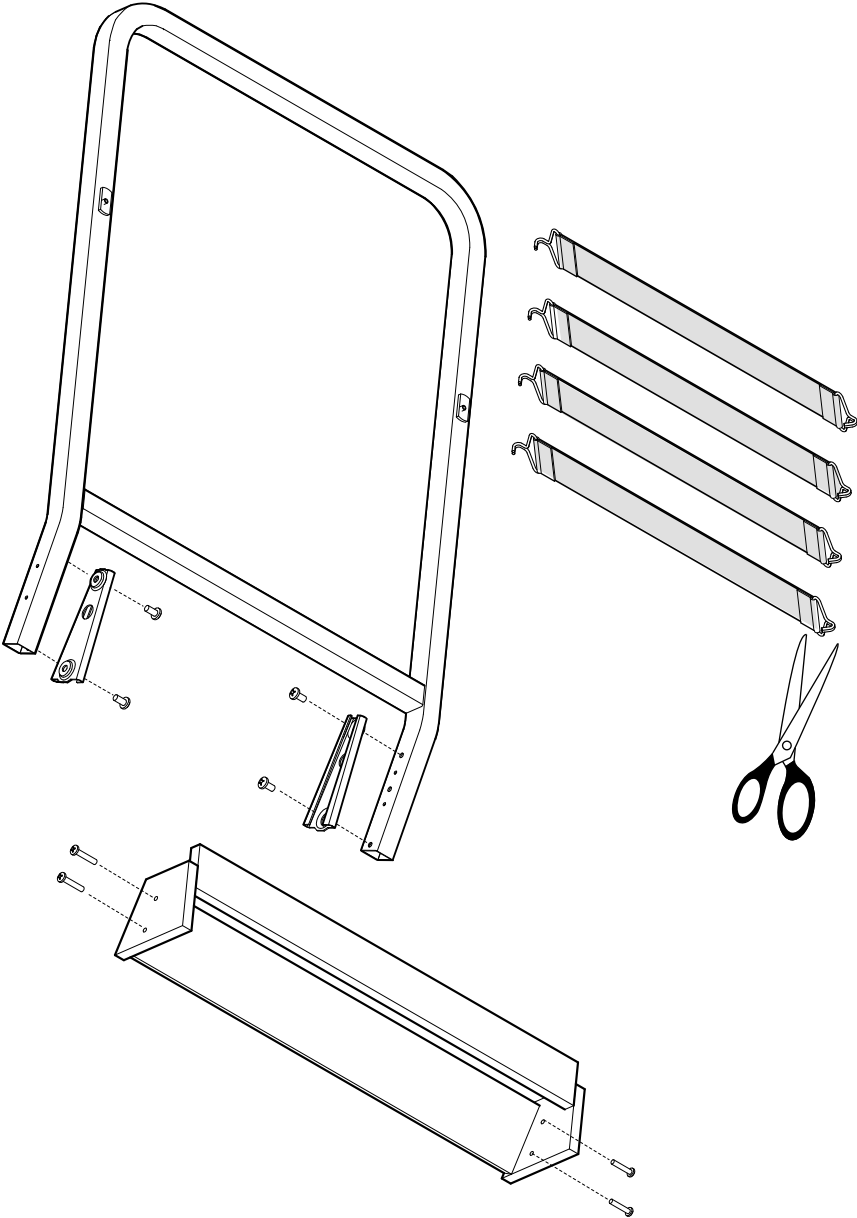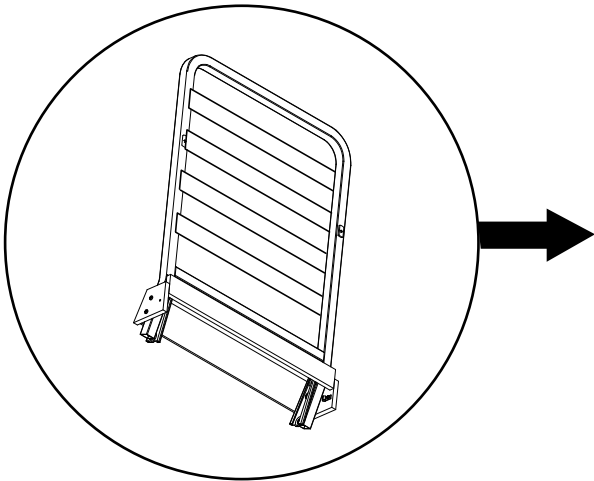

Bins for components removed		
Steel	Fabric	Scrap

Remove Cushions

- 1. Cut off upholstery.
- 2. Disassemble cushions as shown.

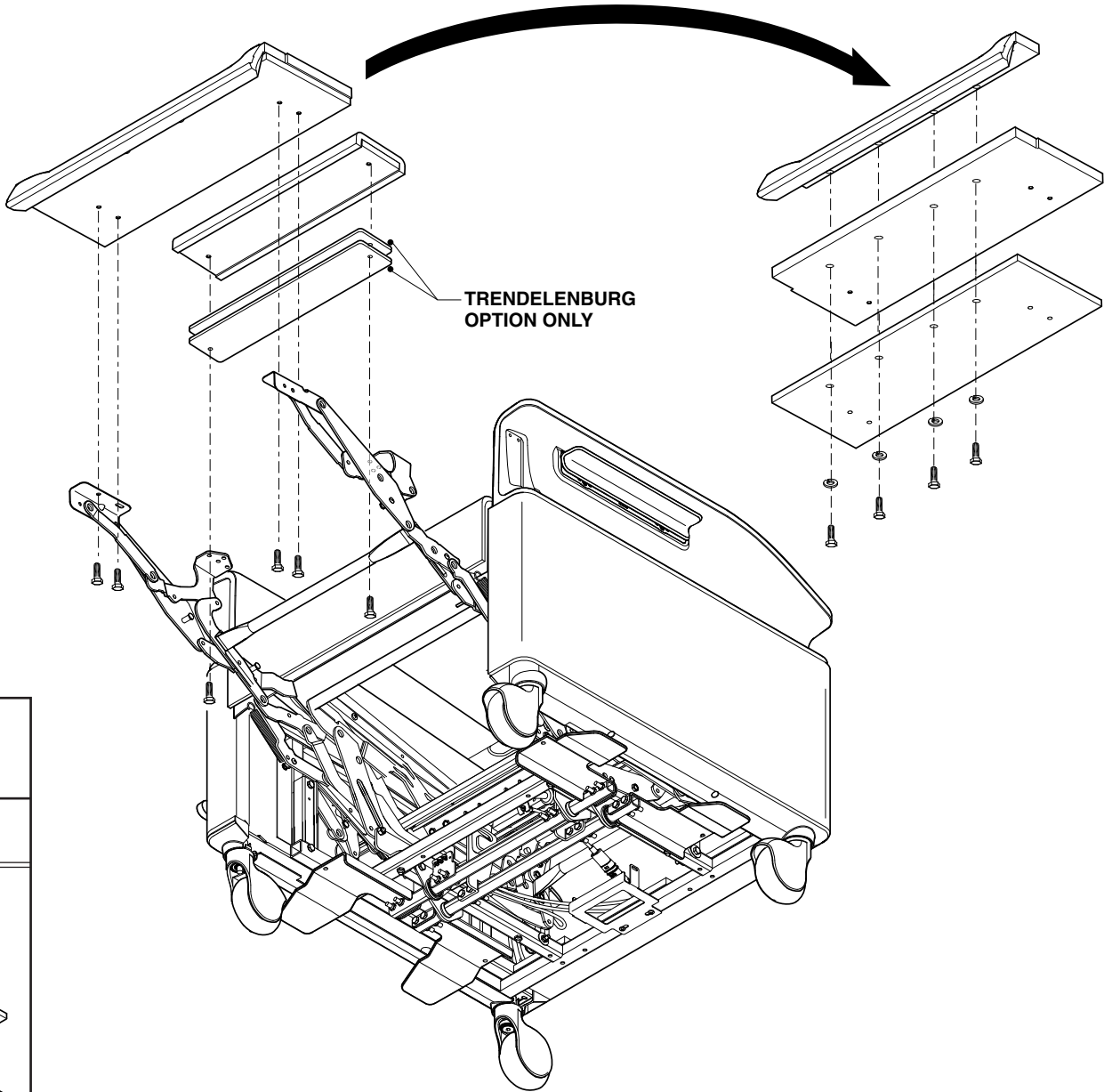


Bins for components removed		
Steel	Upholstery	Scrap

Disconnect Control Cables


Remove Cable Clamps

Remove Handle

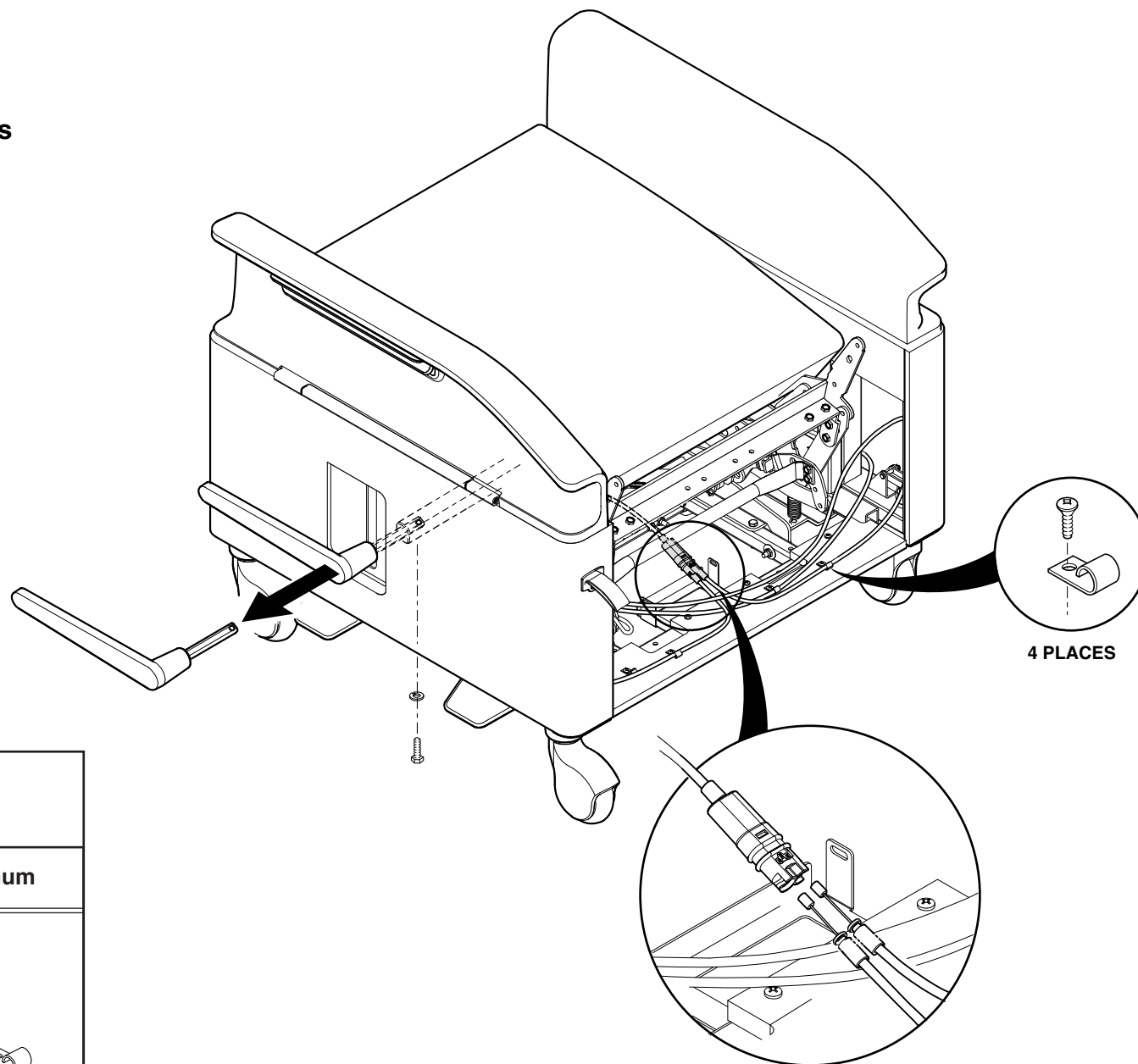


Bins for components removed



Steel

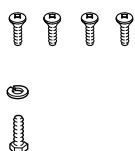

Nylon & Aluminum

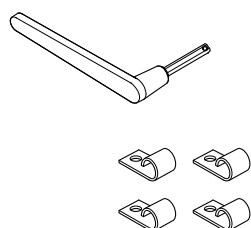


Remove Mechanism

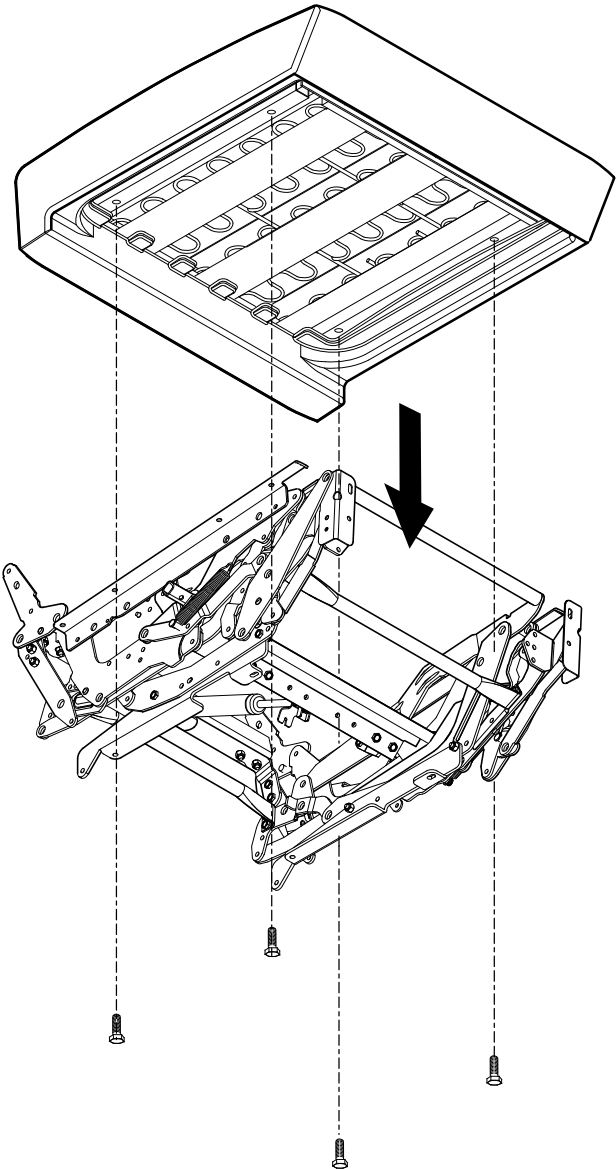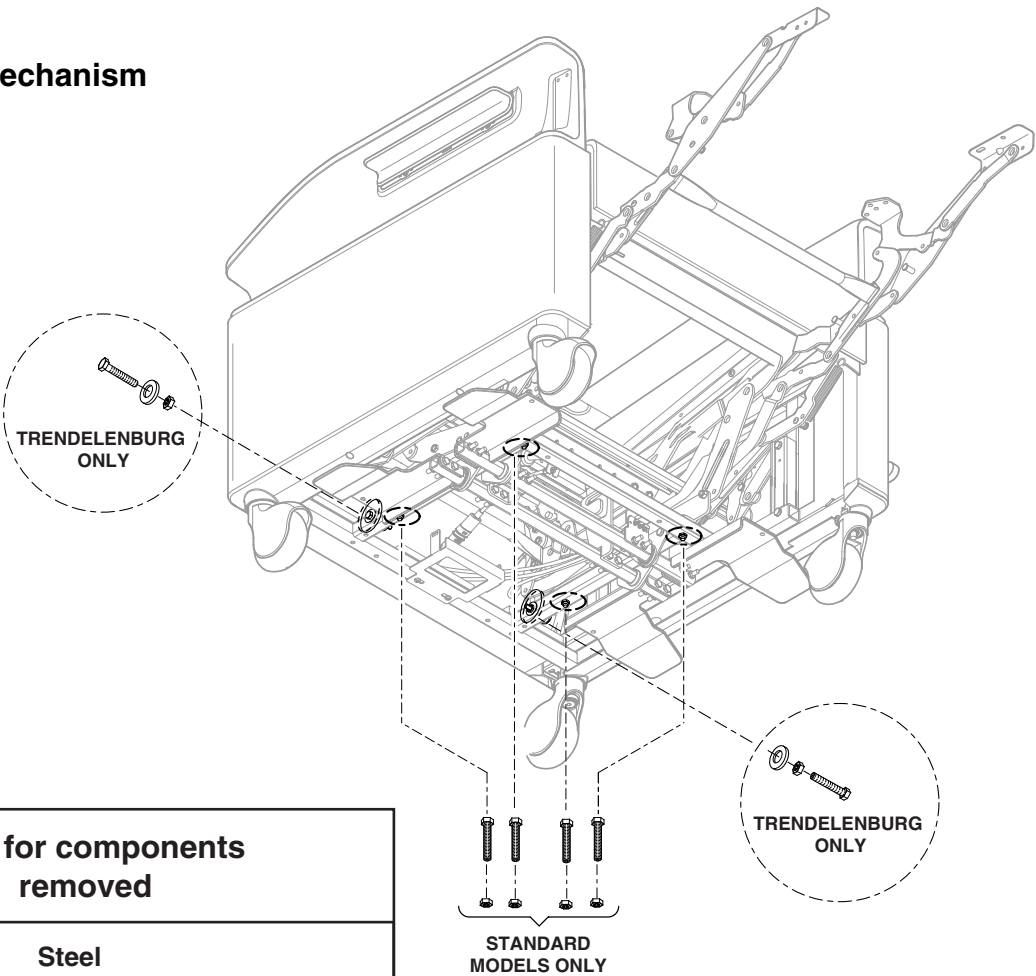


Bins for components removed	
Steel	
	Standard
	Trendelenburg


Cut Upholstery
Disassemble Seat Assembly

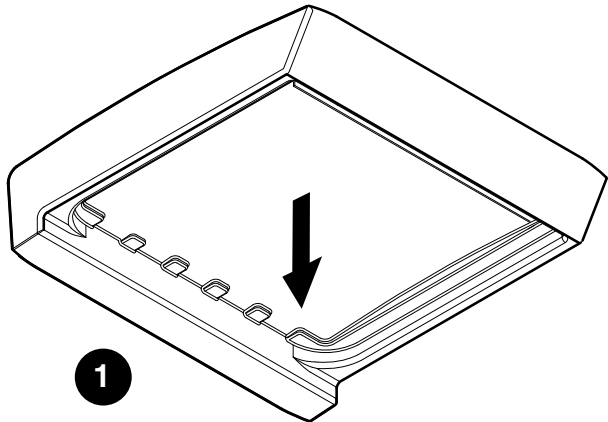


1

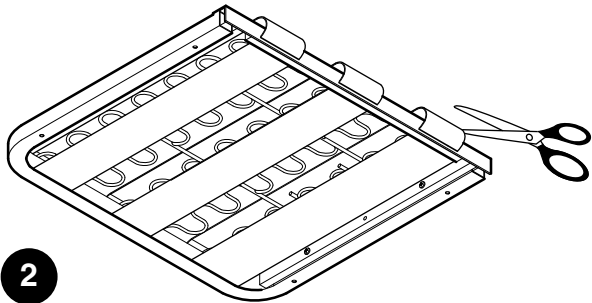





2

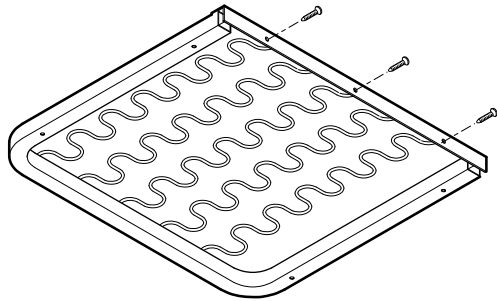



3

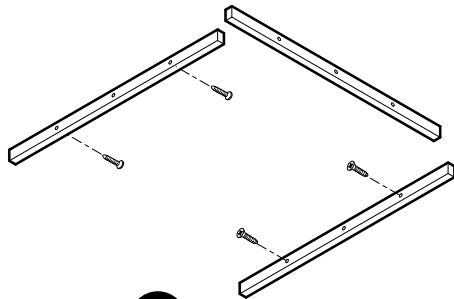

4

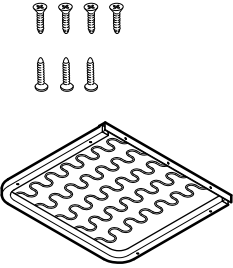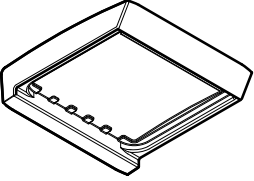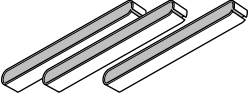
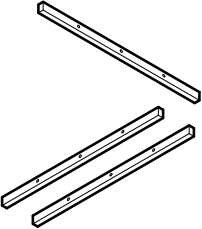
Bins for components removed			
Steel	Foam	Polyester	Scrap
			

**Remove & Disassemble
Fixed Arm**

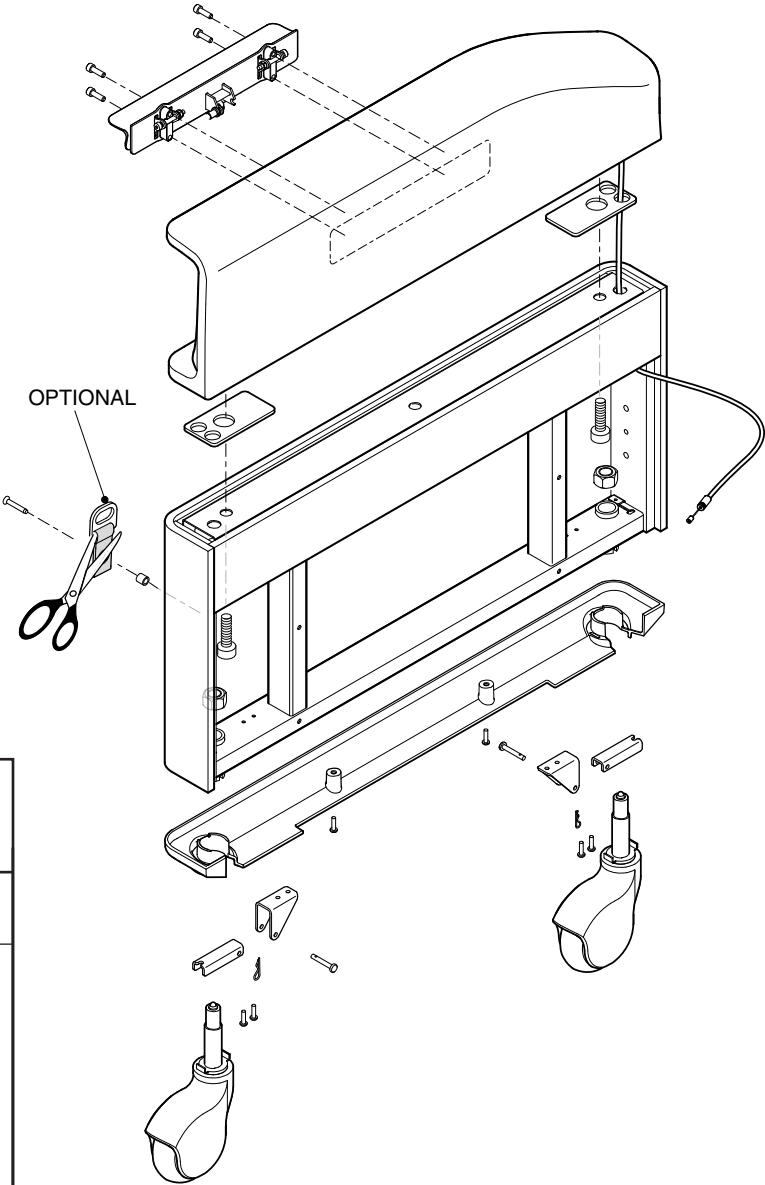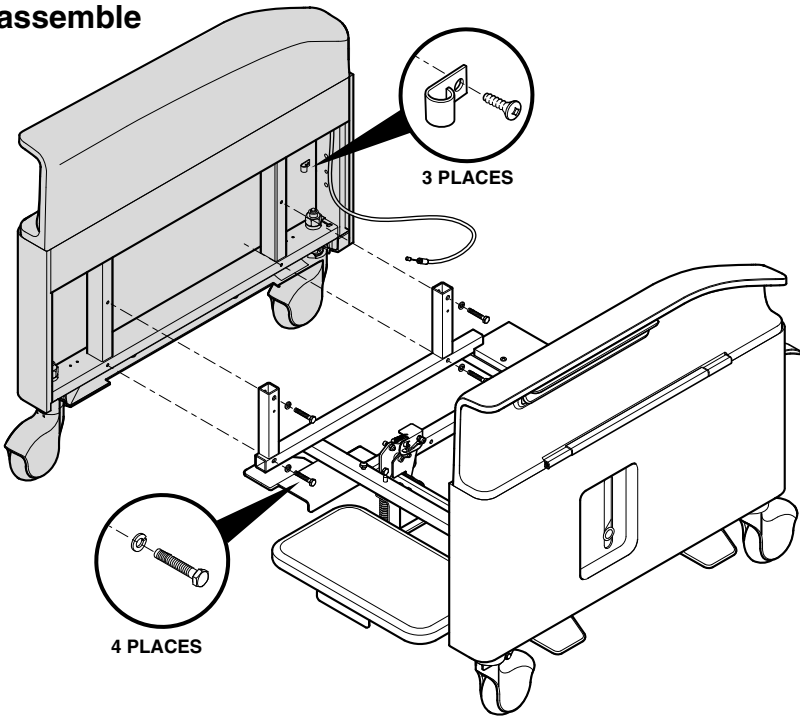


Bins for components removed					
Steel	Aluminum	Polypropylene	Nylon	Bronze	Scrap

**Remove & Disassemble
Folding Arm**

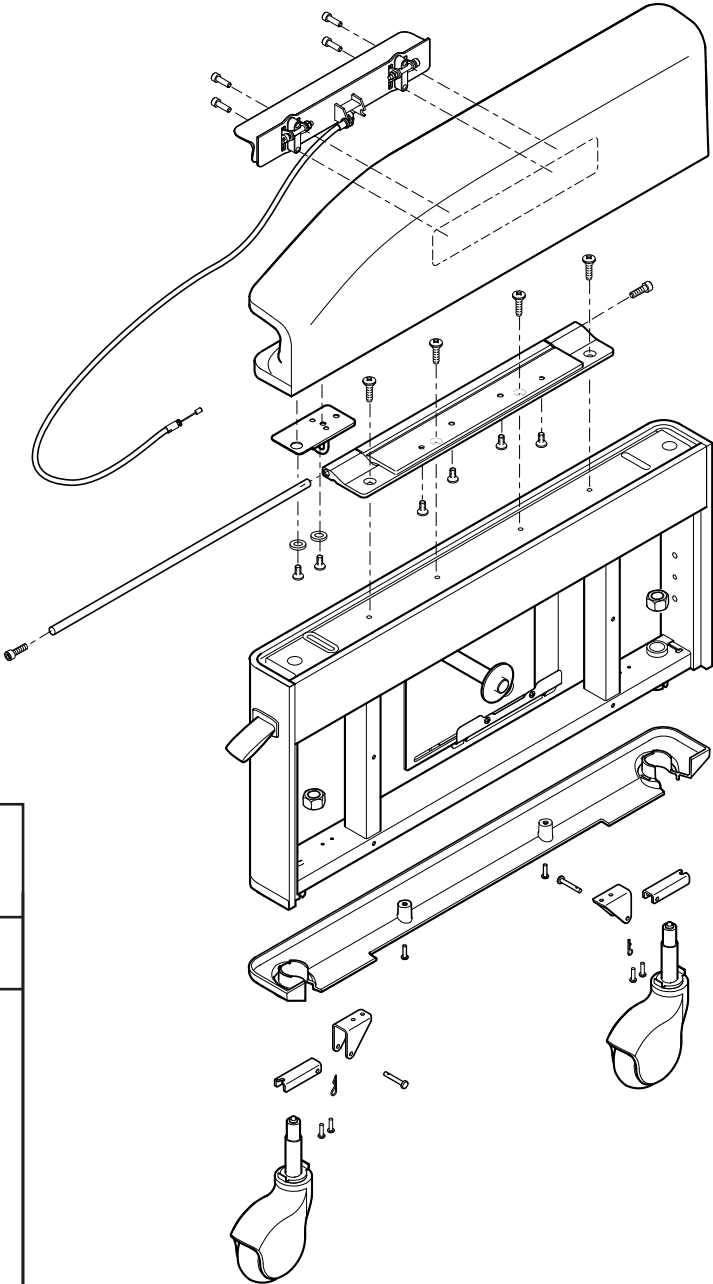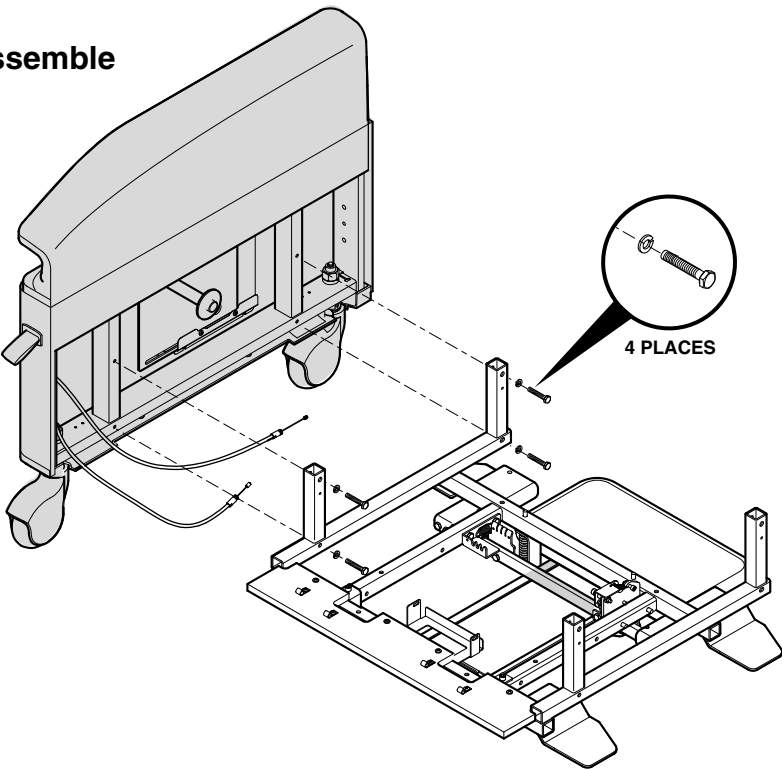


Bins for components removed				
Steel	Aluminum	Nylon	Polypropylene	Scrap



**Disassemble
Lower Arm**

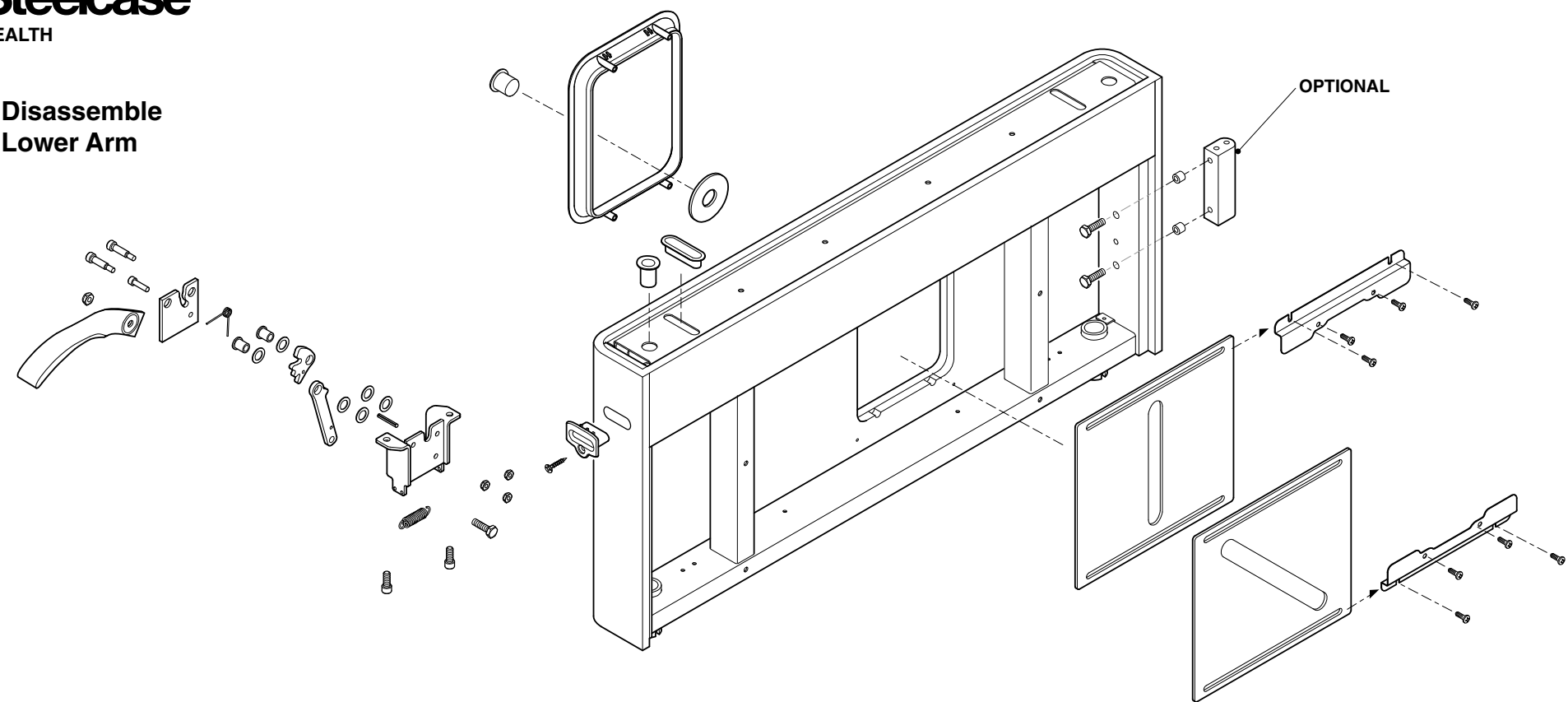


Bins for components removed								
Steel	Polypropylene	HDPE	Zinc	Alum.	Bronze	Acetyl	Nylon	Scrap



Disassemble Frame

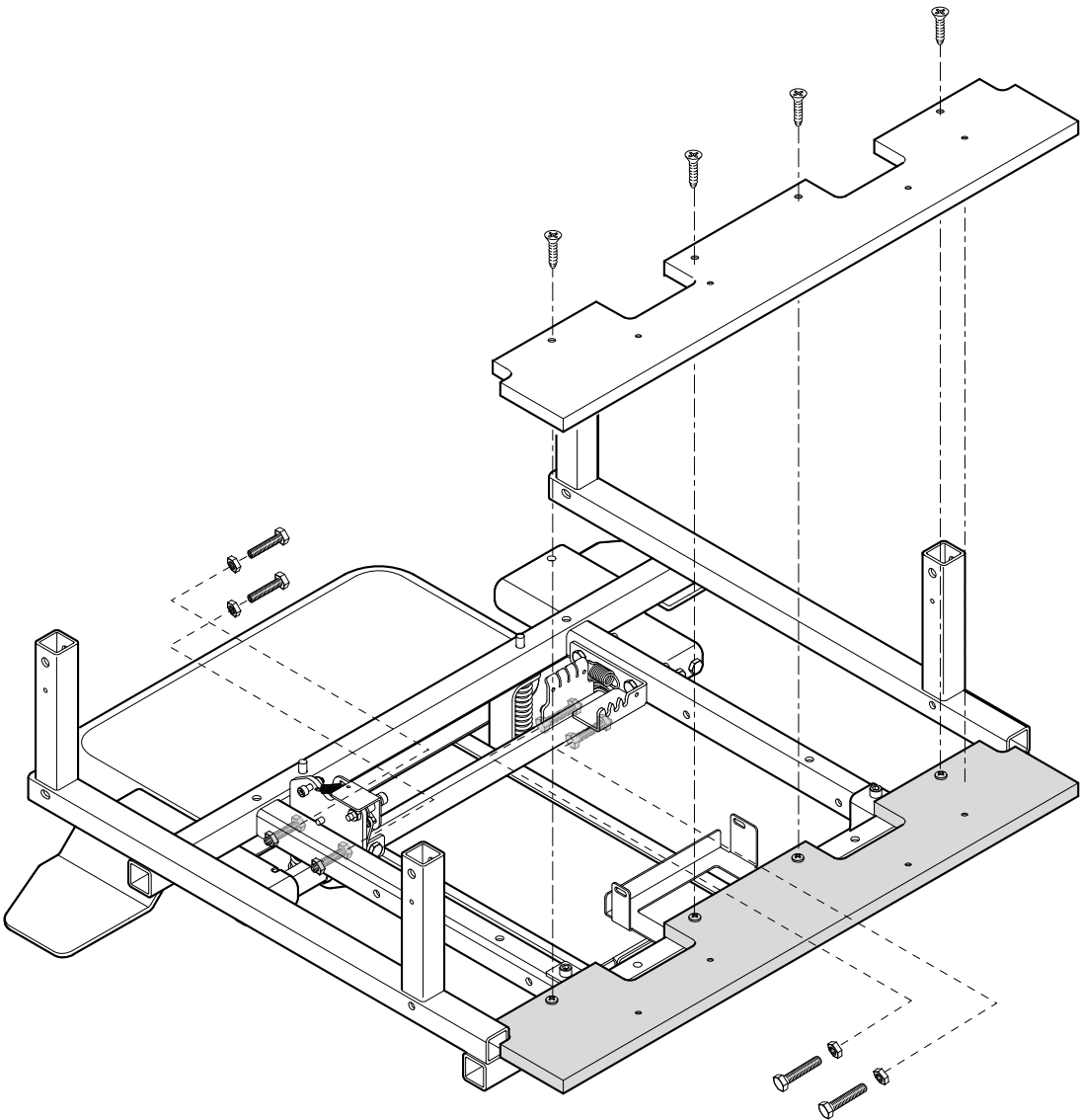

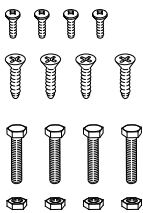
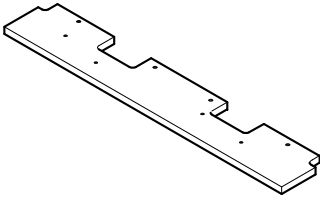
Bins for components removed	
Steel	Scrap
	


**Remove & Disassemble
Optional Footrest**

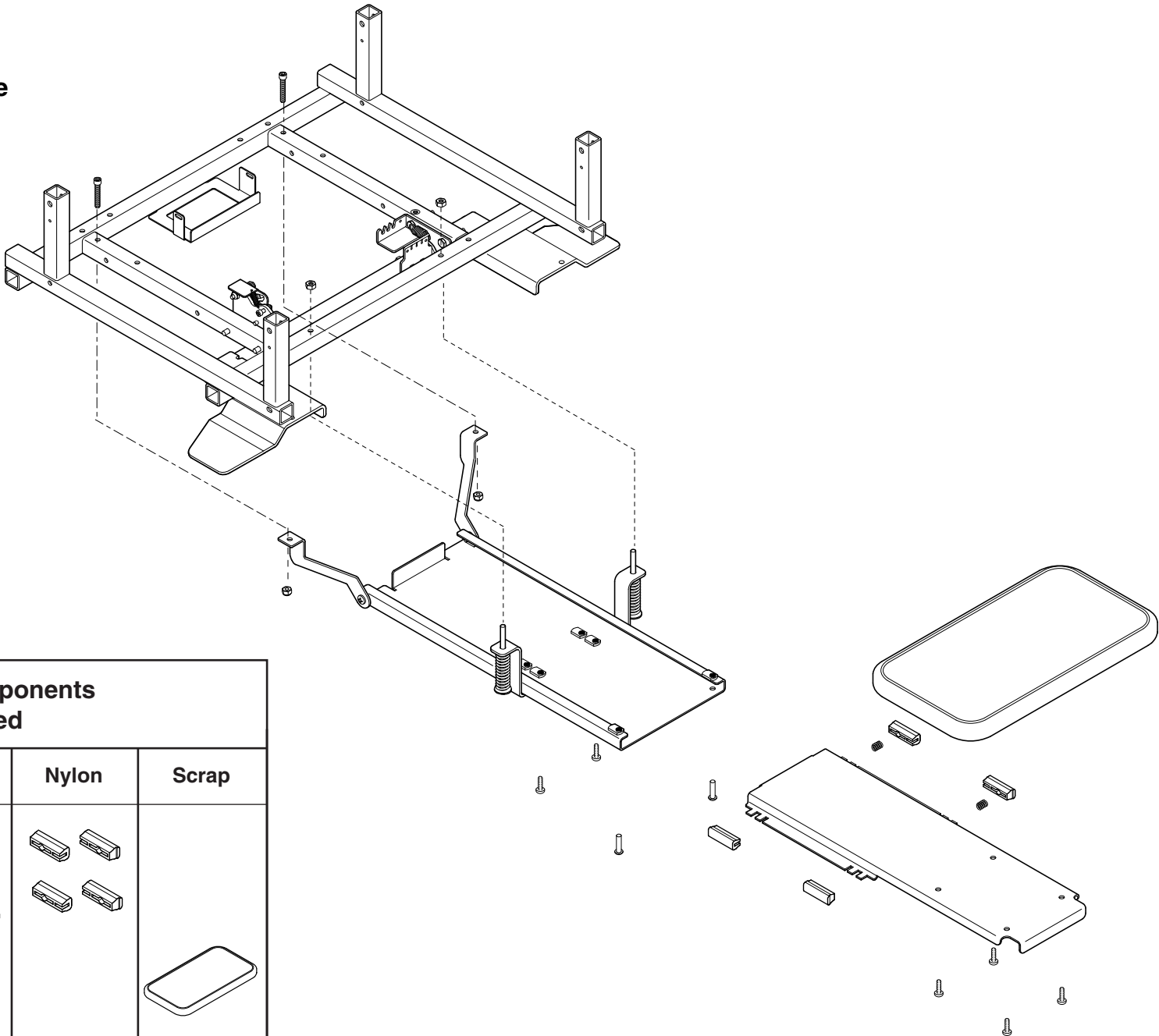

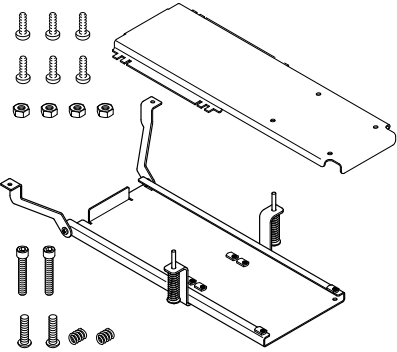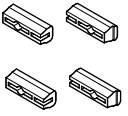
Bins for components removed		
Steel	Nylon	Scrap
		


**Remove & Disassemble
Trendelenburg Option**

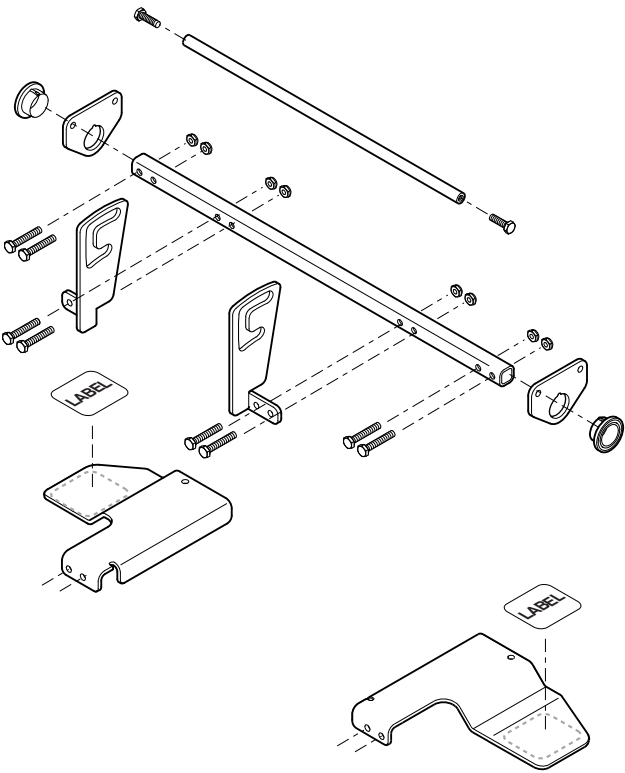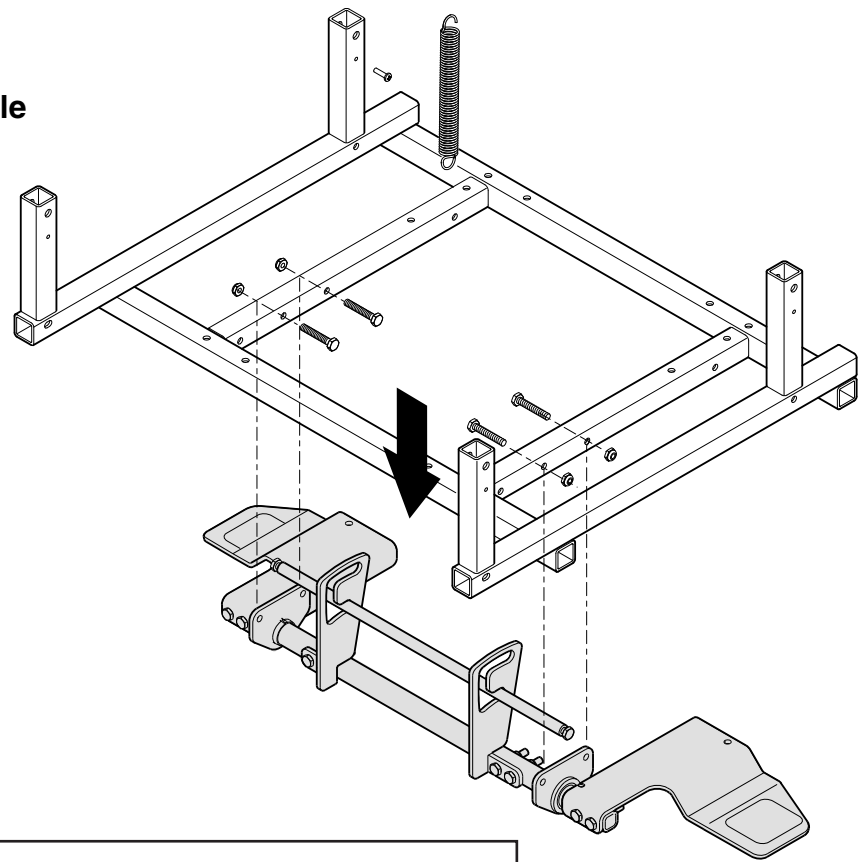

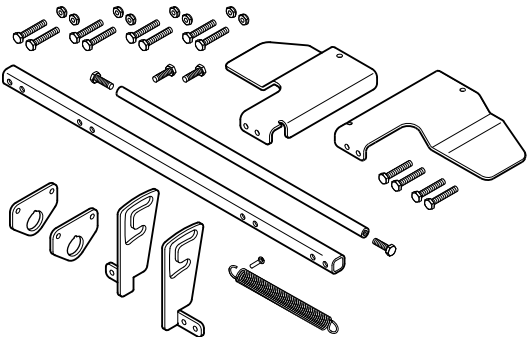
Bins for components removed		
Steel	Nylon	Scrap
		

Remove & Disassemble Brake

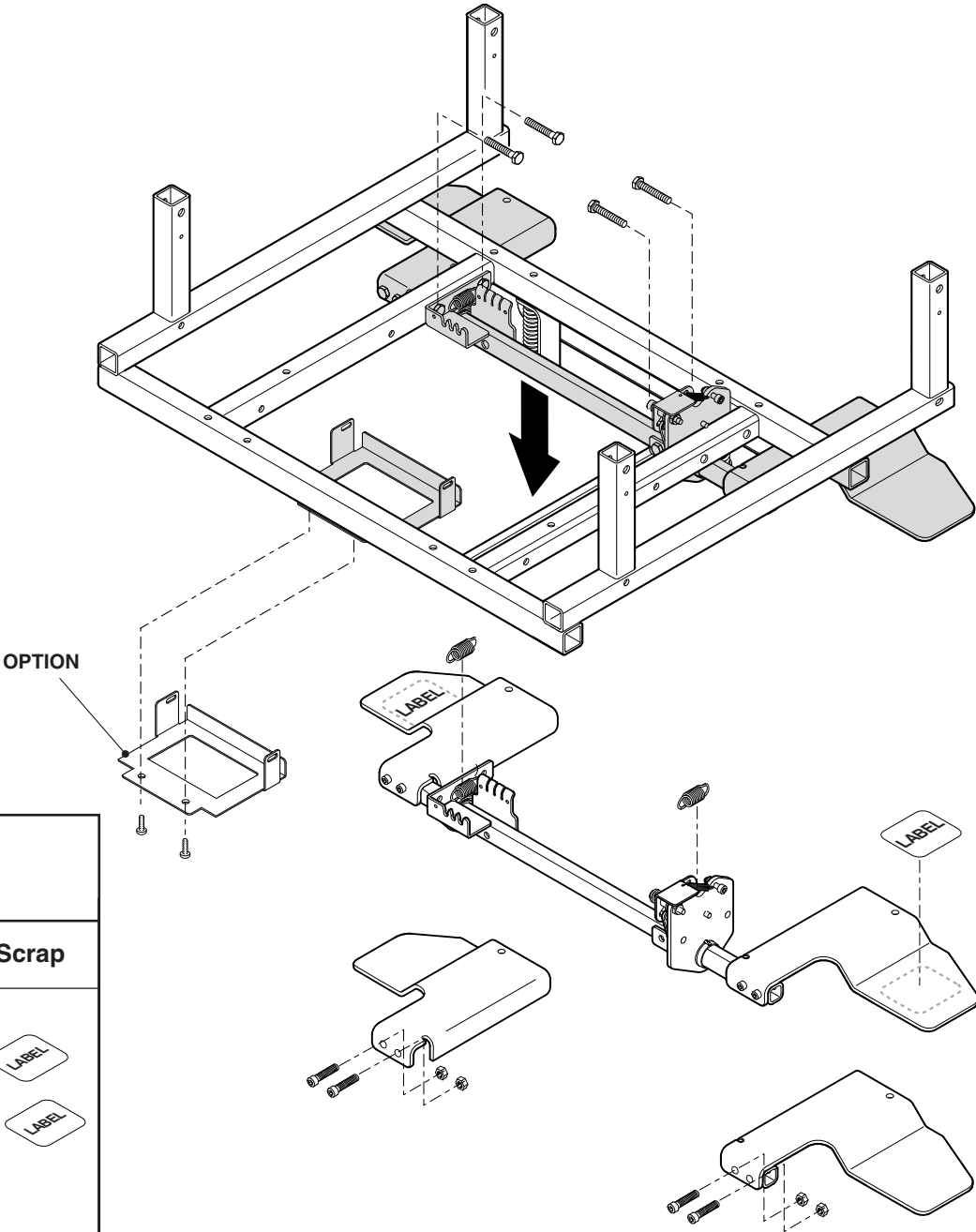

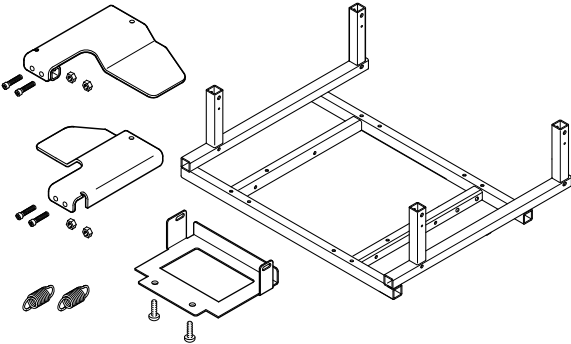
Bins for components removed	
Steel	Scrap
	

Disassemble Brake Bar

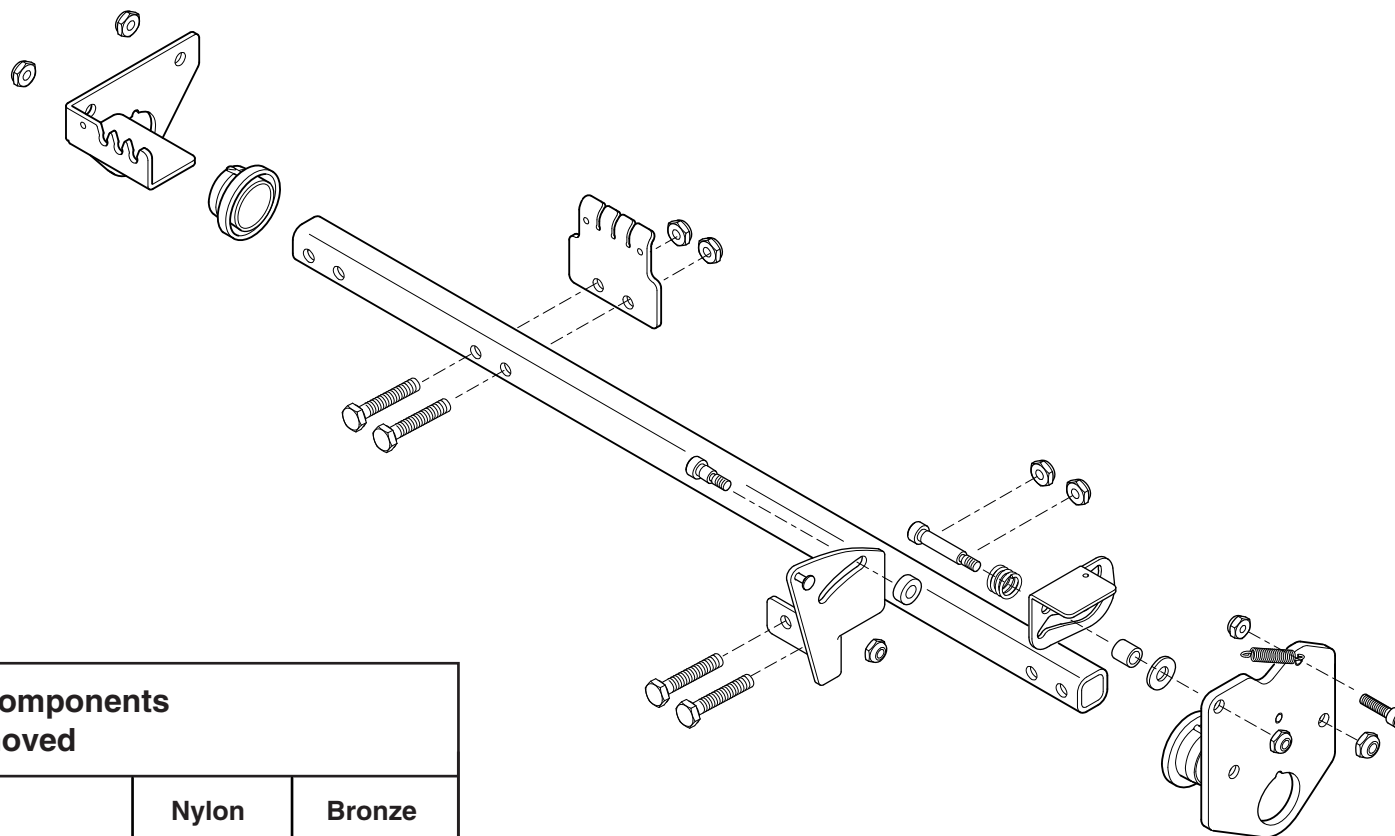

Bins for components removed		
Steel	Nylon	Bronze
A collection of steel components including various bolts, nuts, washers, a long bar, a spring, and several brackets.	Two O-rings and one small circular component.	One small cylindrical component.

**Cut Fabric from Foam
Remove Upholstery,
Foam & Wood**

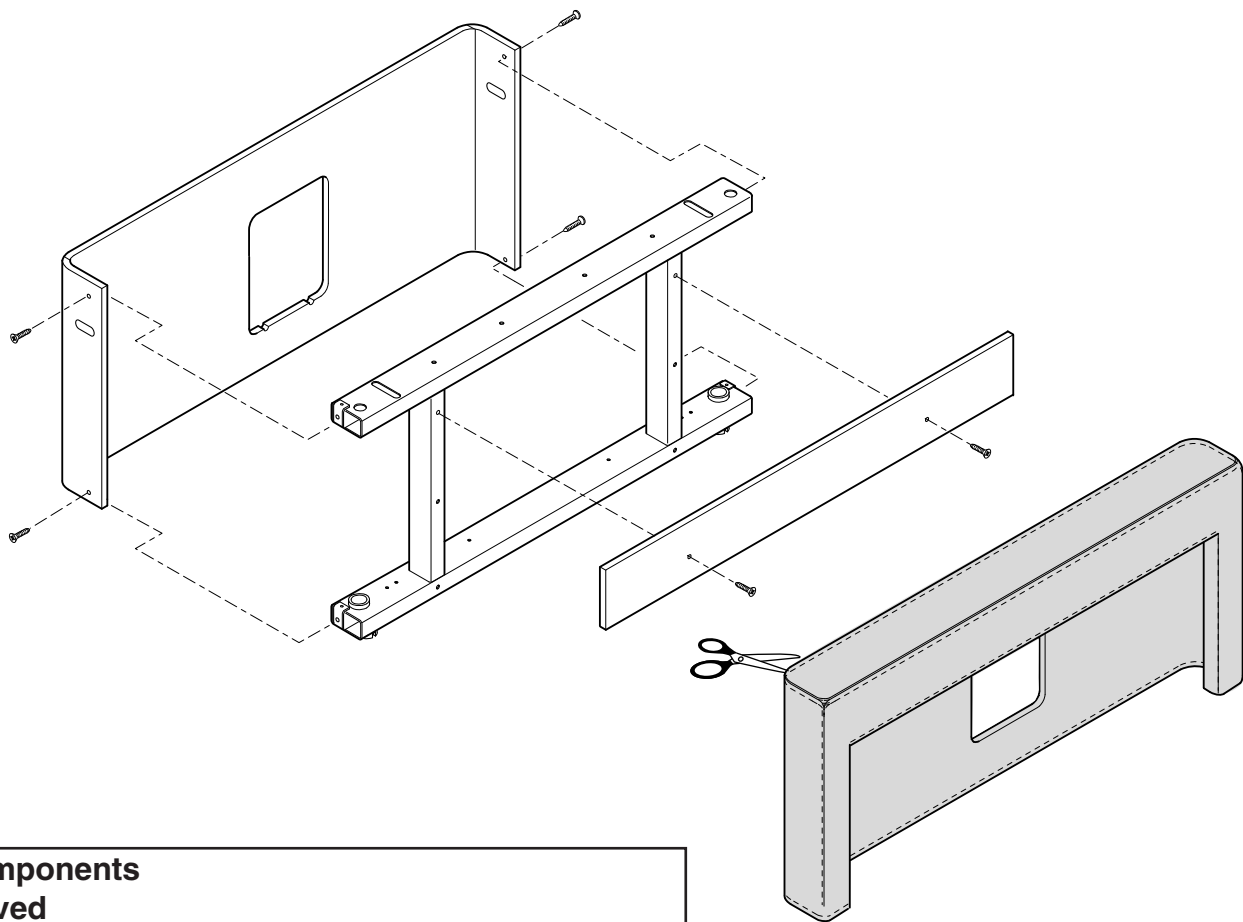

Bins for components removed			
Steel	Foam	Upholstery	Scrap
