


TekTone® Tek-CARE®120 System Layout

(See installation manuals for complete wiring information.)

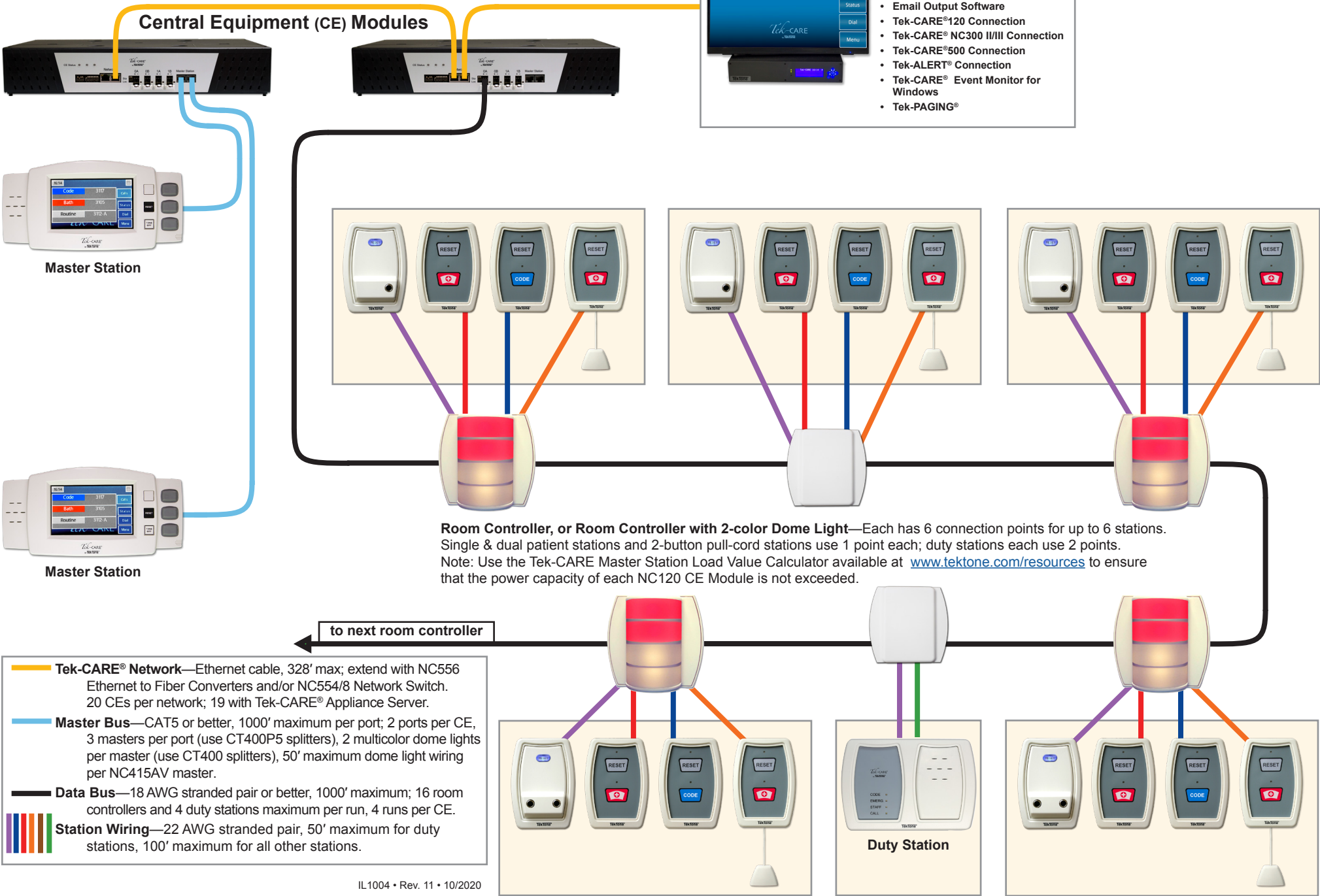

TekTone® Tek-CARE®160 System Layout


(See installation manuals for complete wiring information.)

Central Equipment (CE) Modules

Tek-CARE® Appliance Server Options




- Tek-CARE® Reporting Software
- Tek-CARE TV Display Monitors
- Tek-CARE® Staff App server
- Email Output Software
- Tek-CARE®120 Connection
- Tek-CARE®300 II/III Connection
- Tek-CARE®500 Connection
- Tek-ALERT® Connection
- Tek-CARE® Event Monitor for Windows
- Tek-PAGING®

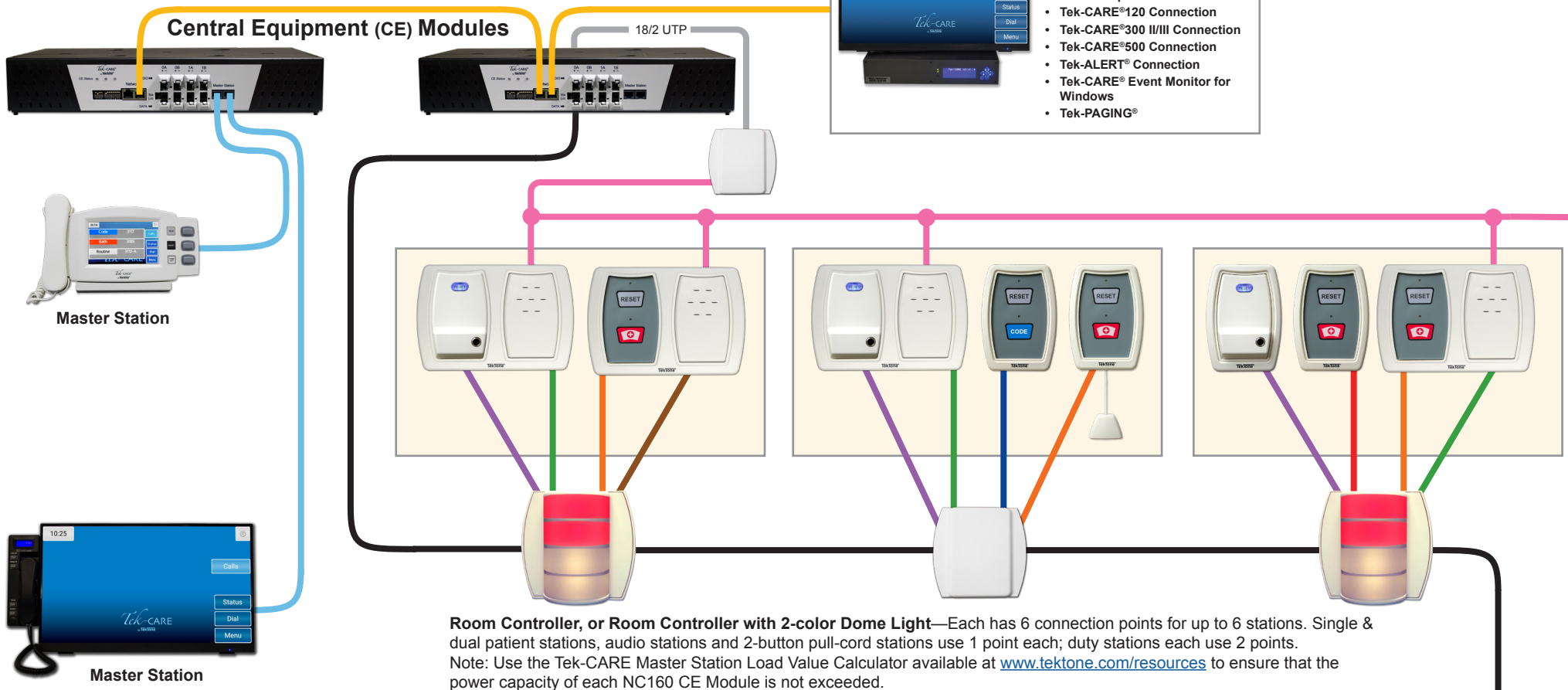

Room Controller, or Room Controller with 2-color Dome Light—Each has 6 connection points for up to 6 stations. Single & dual patient stations, audio stations and 2-button pull-cord stations use 1 point each; duty stations each use 2 points.
 Note: Use the Tek-CARE Master Station Load Value Calculator available at www.tektone.com/resources to ensure that the power capacity of each NC160 CE Module is not exceeded.

- Tek-CARE® Network**—Ethernet cable, 328' max; extend with NC556 Ethernet to Fiber Converters and/or NC554/8 Network Switch. 20 CEs per network; 19 with Tek-CARE® Appliance Server.
- Master Bus**—CAT5 or better, 1000' maximum per port; 2 ports per CE, 3 masters per port (use CT400P5 splitters), 2 multicolor dome lights per master (use CT400 splitters), 50' maximum dome light wiring per NC415AV master.
- Data Bus**—18 AWG stranded pair or better, 1000' maximum; 16 room controllers and 4 duty stations maximum per run, 4 runs per CE.
- Audio Bus**—CAT5e or better, 1000' maximum. 4 runs per CE.
- Station Wiring**—22 AWG stranded pair, 50' maximum for duty stations, 100' maximum for all other stations.
- RS232 Serial Cable**—50' maximum.


IL 1025 R5 02/2022

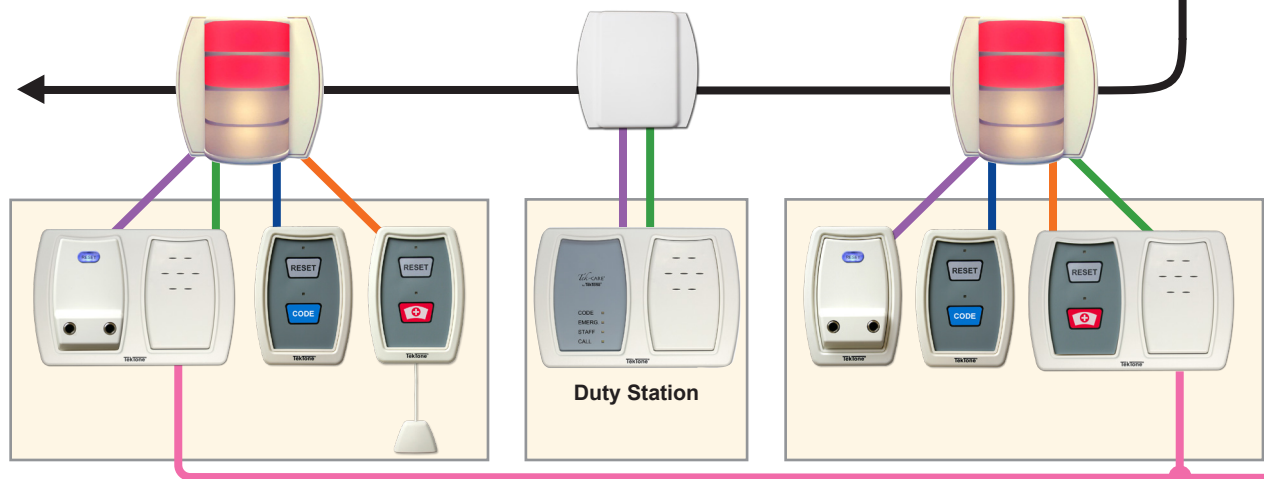

TekTone® Integrated Wired and Wireless Call Systems w/ Domes

IL1091 • Rev. 0 • 02/2022

Tek-CARE120/160
Nurse Call System (Wired)

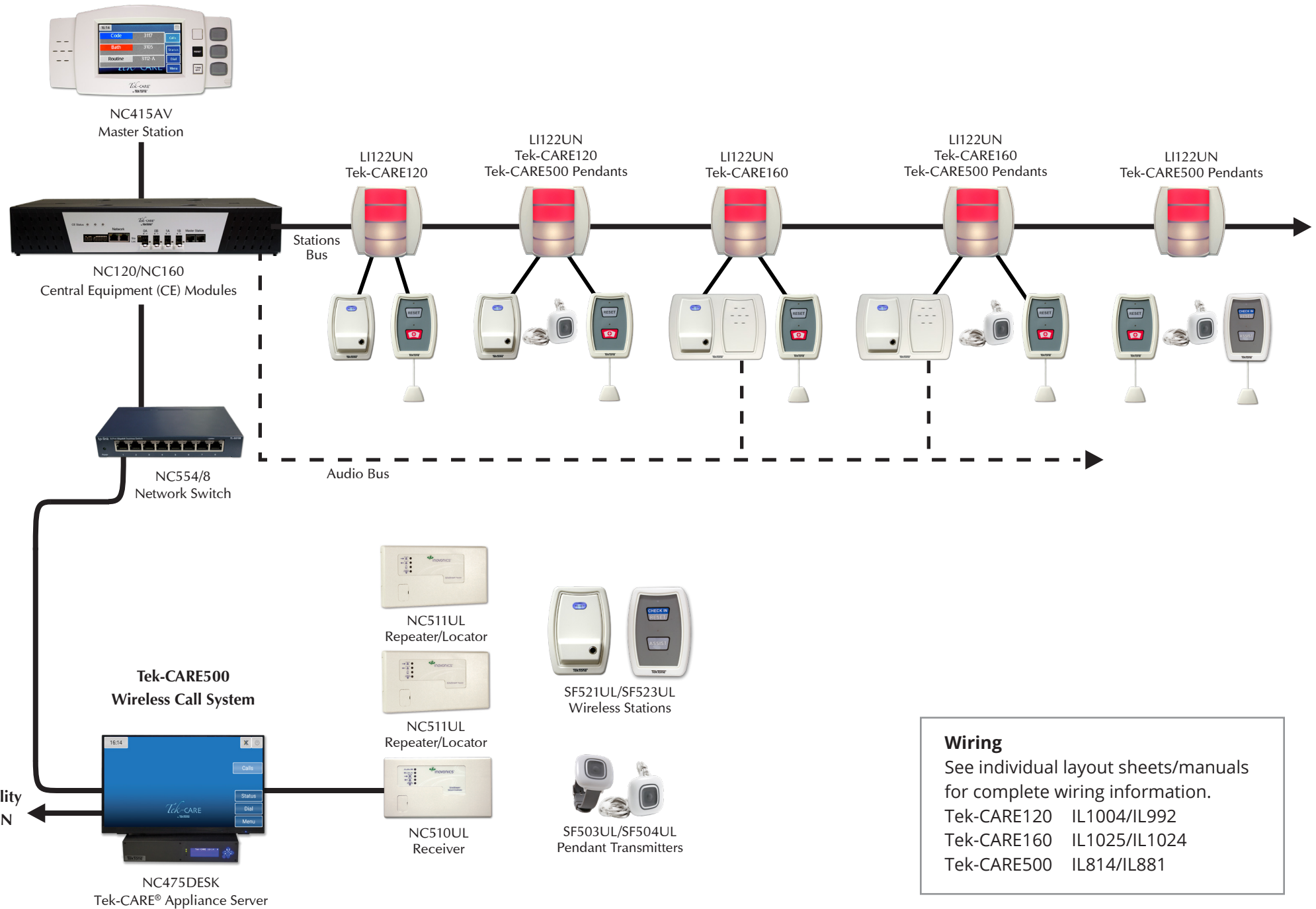

Wiring

See individual layout sheets/manuals for complete wiring information.

Tek-CARE120	IL1004/IL992
Tek-CARE160	IL1025/IL1024
Tek-CARE500	IL814/IL881

TekTone® Tek-CARE® Appliance Server Integration Layout

(See installation manuals for complete wiring information.)

IL919 Tek-CARE Appliance Server color system layout • Rev. 18 • 2/21

Tek-CARE500 Wireless Nurse Call



Tek-CARE400 P5+ Nurse Call



Tek-CARE300 III Nurse Call



Tek-CARE NC110, NC150, NC200 Nurse Call



Tek-CARE160 Nurse Call



Tek-CARE120 Nurse Call



Tek-CARE700 Wander Management



The Tek-CARE Appliance Server can simultaneously connect to any TekTone nurse call system, plus a variety of foreign systems—such as fire alarm and security systems. The Tek-CARE Appliance Server displays the calls, events and messages from all connected systems at once.

NC475 Tek-CARE Appliance Server



Any computer networked with the Tek-CARE Appliance Server can create customizable reports on all activity. Plus, networked PCs can view all system events and messages as they occur via TekTone's Tek-CARE Remote Event Monitor.

Tek-ALERT Integration Manager



Tek-PAGING Radio Pocket Paging



Tek-CARE Remote Event Monitor



Reporting System and Config Tool Live



Media Gateway License




Email Output Software



Staff App iPhone, iPad, iPod Touch, Android



Tek-CARE TV Monitors

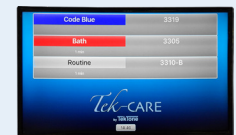

Facility
LAN